

AD-A052 766

HARRY DIAMOND LABS ADELPHI MD  
SKYNET APPLICATIONS SOFTWARE PACKAGE (SASP) OPERATORS' AND USER--ETC(U)  
MAR 78 W J SCOTT, J R HILLER, R PUTTCAMP

F/G 9/2

UNCLASSIFIED

HDL-TR-1843

NL


| OF |  
AD-A052766



END  
DATE  
FILMED  
5-78  
DDC

AD No.             
DDC FILE COPY

AD A 052766

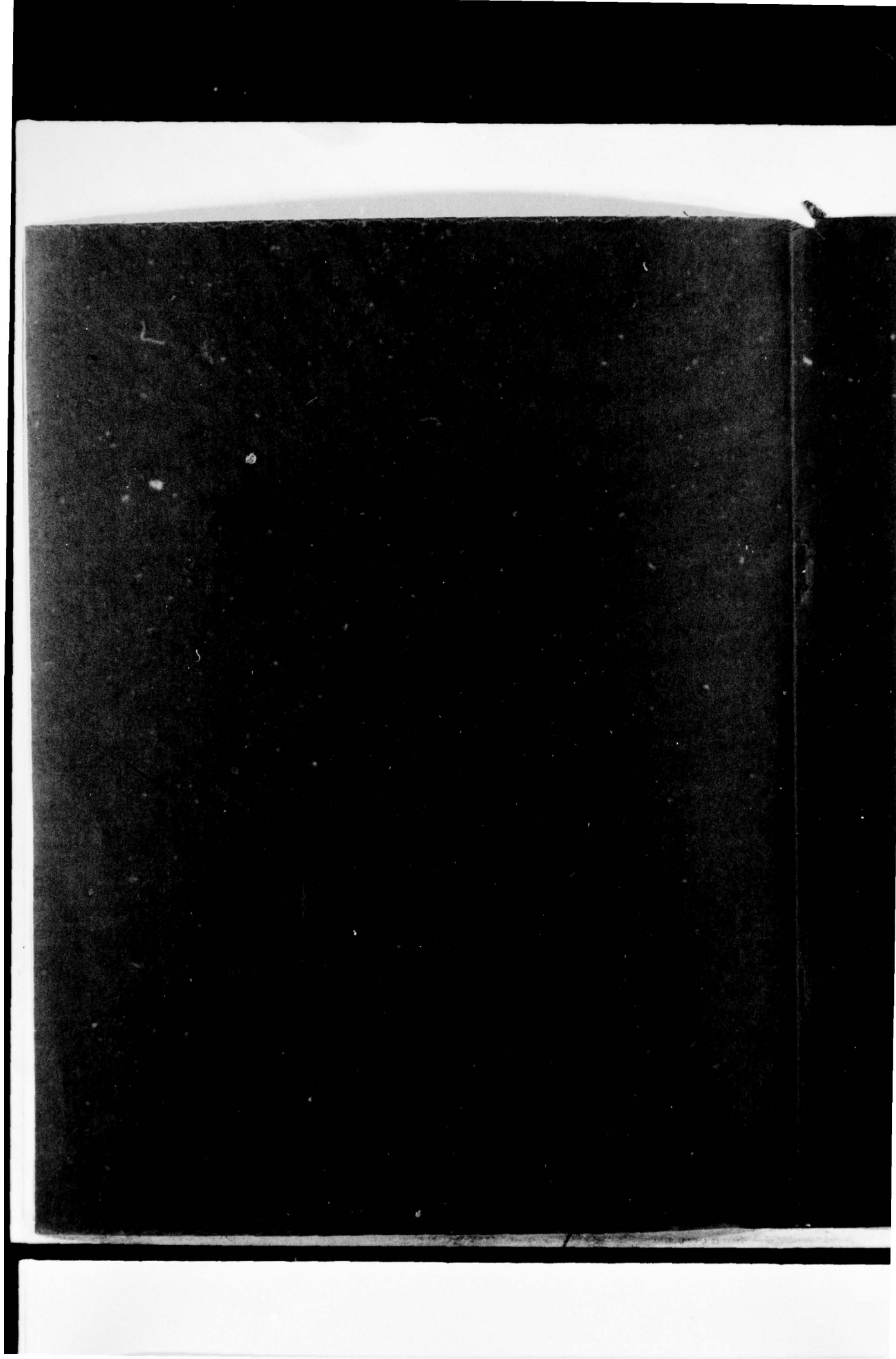

UNCLASSIFIED  
SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

REPORT DOCUMENTATION PAGE		READ INSTRUCTIONS BEFORE COMPLETING FORM
1. REPORT NUMBER 14 HDL-TR-1843	2. GOVT ACCESSION NO.	3. RECIPIENT'S CATALOG NUMBER
4. TITLE (and Subtitle) 6 SKYNET Applications Software Package (SASP) Operators' and Users' Handbook, Version 2	5. TYPE OF REPORT & PERIOD COVERED 9 Technical Report.	
7. AUTHOR(s) 10 Walter J. Scott, John R. Hiller, Robert Puttcamp	6. PERFORMING ORG. REPORT NUMBER	
8. CONTRACT OR GRANT NUMBER(s)	10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS Program Ele: 6.27.04.H	
9. PERFORMING ORGANIZATION NAME AND ADDRESS Harry Diamond Laboratories 2800 Powder Mill Road Adelphi, MD 20783	11. CONTROLLING OFFICE NAME AND ADDRESS Director Defense Nuclear Agency Washington, DC 20305 11	12. REPORT DATE Mar 1978
14. MONITORING AGENCY NAME & ADDRESS (if different from Controlling Office) 16 R99QAXE 17 E502	13. NUMBER OF PAGES 93 12 89p	15. SECURITY CLASS. (of this report) Unclassified
15a. DECLASSIFICATION/DOWNGRADING SCHEDULE		
16. DISTRIBUTION STATEMENT (for this Report) Approved for public release; distribution unlimited.		
17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report)		
18. SUPPLEMENTARY NOTES HDL Project: 284724 DRCMS Code: 36AA.7100.62704 This research was sponsored by the Defense Nuclear Agency under Subtask R99QAXEEB502, Systems Generated EMP, Work Unit 06, SKYNET.		
19. KEY WORDS (Continue on reverse side if necessary and identify by block number) SGEMP SKYNET Data acquisition software Minicomputer applications software		
20. ABSTRACT (Continue on reverse side if necessary and identify by block number) Data acquisition and reduction on a large-scale systems generated electromagnetic pulse (SGEMP) test program are traditionally costly and tedious. The most widely used methods of acquiring and recording the system response waveforms involve the use of oscilloscopes and cameras, the waveforms being captured on Polaroid film. This method requires that the analog waveforms		

DD FORM 1 JAN 73 1473 EDITION OF 1 NOV 65 IS OBSOLETE

UNCLASSIFIED

1 SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

163 050 - JOB


UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE(When Data Entered)

(oscillograms) be digitized, usually by hand, into a computer-compatible format prior to applying automated data analysis techniques. The amount of manual data handling and the need for many skilled technicians make this an error-sensitive and costly procedure.

This report presents a new procedure for capturing and processing complex SGEMP transients. The method is incorporated into an SGEMP Transportable Automated Recording System (STARS). It reduces manual data handling significantly (a major source of errors); outputs data as computer compatible magnetic tapes, annotated hard copies, or both; and lessens the manpower requirements for data acquisition and processing. The software implementation, SKYNET Applications Software Package (SASP), is discussed here to provide STARS operators and users with the detailed information needed to effectively utilize the system. A companion volume, the SKYNET Applications Software Package (SASP), Version 2, Programmers' Handbook, contains detailed summaries of all SASP source codes.

This report is a revision of HDL-TR-1817, SKYNET Applications Software Package (SASP) Operator's and User's Handbook, by Walter J. Scott and John R. Hiller (August 1977). This report reflects major changes and additions that have been incorporated into SASP.

ACCESSION FOR	
RTIS	Write Section <input checked="" type="checkbox"/>
ODS	Diff Section <input type="checkbox"/>
UNANNOUNCED	<input type="checkbox"/>
JUSTIFICATION	
BY	
DISTRIBUTION/AVAILABILITY CODES	
Dist.	AVAIL. CODE OR SPECIAL
	

UNCLASSIFIED

## CONTENTS

	<u>Page</u>
1. INTRODUCTION . . . . .	7
1.1 Background . . . . .	7
1.2 Overview . . . . .	7
2. DATA FILES . . . . .	11
2.1 Introduction . . . . .	11
2.2 File Structures . . . . .	13
3. DATA SYSTEM OPERATIONS . . . . .	22
3.1 Introduction . . . . .	22
3.2 Data System Generation . . . . .	23
3.3 Terminal Input Conventions . . . . .	33
3.4 Data-Taking Sequence . . . . .	33
3.4.1 TESTSEQ Operations Control File .. . . .	33
3.4.2 Manual Data Acquisition Procedures . . . . .	42
3.4.3 Disc Preparation . . . . .	42
4. DATA PROCESSING . . . . .	42
4.1 Normalization . . . . .	42
4.2 Data Evaluation and Characterization . . . . .	45
4.2.1 Baseline Data . . . . .	46
4.2.2 Calibration Data . . . . .	46
4.2.3 Raw Data . . . . .	47
4.2.4 Normalized Data . . . . .	48
4.3 Typical Production Run . . . . .	49
4.4 Storage and Retrieval . . . . .	52
5. BRIEF MODULE DESCRIPTIONS . . . . .	53
ABBREVIATIONS . . . . .	61
DISTRIBUTION . . . . .	91

## APPENDICES

	<u>Page</u>
A.--DIGITIZER AND CHANNEL ERROR CODES . . . . .	63
B.--DIGITIZER COMMANDS AND STATUS . . . . .	69
C.--TESTSEQ OPERATIONS CONTROL FILE . . . . .	73
D.--SYMDEF MACRO ENTRIES . . . . .	79
E.--LIBRARIES: OBJECT AND MACRO . . . . .	89

## FIGURES

1 Systems generated electromagnetic pulse Transportable Automated Recording System (STARS) . . . . .	10
2 Calibration file structures . . . . .	17-18
3 System configuration (SYSCON.DAT) file structure . . . . .	19-20
4 Shot file structure . . . . .	21-22
5 Generalized data channel . . . . .	25
6 Delay considerations . . . . .	27
7 Probe log . . . . .	29
8 Transmitter log . . . . .	30
9 Receiver log . . . . .	31
10 Digitizer log . . . . .	32
11 Flowchart for TESTSEQ . . . . .	35-41
12 Raw data, sample 1 . . . . .	47
13 Normalized data . . . . .	48
14 System configuration printout . . . . .	49
15 Raw data, sample 2 . . . . .	50
16 Scaled data . . . . .	51
17 Probe file listing . . . . .	51

## TABLES

	<u>Page</u>
I    Data Files . . . . .	12
II   Calibration File Entries . . . . .	14
III   System Configuration (SYSCON.DAT) File Entries . . . . .	15
IV   Main Data (SHOT.DAT) File Entries . . . . .	16

## 1. INTRODUCTION

### 1.1 Background

In a systems generated electromagnetic pulse (SGEMP) test program, time-history representations of electrical responses at various critical nodes within a test object are provided to analysts for evaluation. Past methods of acquiring and recording these response waveforms involved the use of oscilloscopes and cameras, the waveforms being captured on Polaroid film. Analysis was frequently conducted by use of a digital computer. This analysis required that the analog waveforms be digitized into a computer-compatible digital format. Several types of errors, as well as the time delay in processing the data, are inherent in the manual digitization of oscillographs. For example, as the digitizer operator becomes fatigued (after digitizing several oscillograms), he is likely to move the cursor in such a way that extraneous high frequencies are introduced into the waveform being digitized. The amount of error introduced is directly relatable to the oscillograph trace width (bloom) and the trace derivative<sup>1</sup> (steepness of slope).

This report presents a new method for capturing and processing complex SGEMP transients. The method is incorporated into an SGEMP Transportable Automated Recording System (STARS). It reduces manual data handling significantly (a major source of errors); outputs data as computer-compatible magnetic tapes, annotated hard copies, or both; and lessens the manpower requirements for data acquisition and processing personnel. Thus, the STARS enhances SGEMP data accuracy and reliability through elimination of hand digitization and through better signal recording resolution, while substantially reducing test operational costs. This report is a direct replacement for HDL-TR-1817, SKYNET Applications Software Package (SASP) Operator's and User's Handbook by Walter J. Scott and John R. Hiller (August 1977). It reflects major changes and additions which have been incorporated into SASP.

### 1.2 Overview

The content of this handbook is directed towards STARS operators and data users, the text being chiefly concerned with the description of the SKYNET Applications Software Package (SASP). A brief synopsis of each section follows.

Section 2 is a complete, detailed presentation of system calibration and data files. A thorough understanding of file

---

<sup>1</sup>Data Reduction Planning and Coordination, EG&G, Albuquerque, NM (1972).

structures--their content and usage--is basic to understanding subsequent sections. Through familiarization with the various file entries, the user can gain an appreciation for what information he must provide the system before testing can take place, as well as for the structure of the data bank.

In section 3, user input requirements are described in detail. Sections 3.2 and 3.3 provide the operator with the knowledge needed to run the STARS console terminal and generate or modify the system files. The main operation and control program, indirect command file TESTSEQ, is presented in section 3.4.1. Modular in construction, SASP is composed of well-defined programs that function independently of one another; however, the order in which they are executed is important. Section 3.4.2 discusses some of the do's and don'ts for running single programs in the manual mode.

The raw data acquired from the transient digitizer is not in a form compatible for processing by digital computer. Usually, there is an overabundance of data describing the waveform. Sometimes, there is a lack of data describing certain portions of the waveform. It is the task of the normalization procedure (sect. 4.1) to process these raw data, with their excesses and deficiencies, and to obtain from the raw data a more usable presentation of the waveform.

The STARS operator is tasked with reviewing the data for quality and completeness. He performs data evaluation (sect. 4.2) independently of the data's context to the experiment. Rather, data are examined to identify hardware relatable problems, e.g., data exceeding the digitizer writing rate. Section 4.3 presents a typical production run including user output plots and tables. Section 4 concludes with a discussion of data storage and retrieval.

Section 5 summarizes each SASP module: the indirect command file, program, and subroutines. These summaries, when used with the flowcharts (sect. 3.4.1) and the discussion that follows, provide an overview of the data taking. A more detailed explication of the modules may be found in the SASP Programmers' Handbook.<sup>2</sup>

The SASP software was developed for two reasons:

a. To control the Harry Diamond Laboratories optical system, the calorimeter data acquisition channels, and the Tektronix R7912 transient digitizers.

---

<sup>2</sup>Robert Puttcamp, *SKYNET Applications Software Package (SASP) Programmers' Handbook, Version 2, Harry Diamond Laboratories TR-1844 (February 1978)*.

b. To record, process, store, and reproduce SGEMP test data.

To appreciate many of the software tasks and to obtain an overview of the entire STARS operation, a brief presentation of the data-taking sequence is now developed.

Figure 1 is a STARS system block diagram. Prior to any data taking, the STARS operator must generate the system files by typing in the user-supplied data for the probes, sensors, devices, cables, etc., to be used in the experiment. This information is input at the main console graphic display Tektronix R4010 cathode ray tube terminal and stored on disc unit No. 0 (DK0:), a Digital Equipment Corp. (DEC) RK11 disc. This is designated the "system" disc and contains the system calibration and configuration files, the DEC RSX-11M operating system and the SASP program files. Disc unit No. 1 contains the shot data files, work files used by the Versatec plotting package and user terminal storage.

Once the data channels are completely defined, they are "known" to the system. Any calculations on the data requiring probe sensitivity factors, time base settings, or any other parameter can be made by retrieving the entry from the appropriate system file. The process of data acquisition can now proceed.

The first step in acquiring data is to adjust the knob settings on the digitizers. The Tektronix R632 television (TV) monitors present a visual display to aid the operator in the adjustment of the vertical and horizontal baseline positions, as well as the intensity settings. The operator can adjust them by running the program KNOB and making the required adjustments. An operations control file, TESTSEQ, itself composed of individual command lines, allows interactive control of the entire data-taking sequence. Then program BASE is run to acquire background data from each digitizer. These data are stored on DK1: in the file named MAIN.DAT. Program BASE may be repeated if errors have been flagged. (Error flags are discussed in app A.) Otherwise, the operator continues by acquiring calibration signals from the calibration generators located inside the optical transmitters. Program CAL accomplishes the entire calibration sequence, including automatic turn on of the optical equipment and triggering of the digitizers. Program CAL data are then plotted on the terminal for evaluation by the STARS operator. Again, the flexible TESTSEQ indirect command file permits looping back to reacquire these calibration waveforms if desired. The next program, PREDAT, prepares the system for the actual data shot.

Program PREDAT sends commands to all digitizers and places them in the proper mode for transient recording. (App B discusses these commands in detail.) The acquisition and storage of raw digitizer data and the acquisition, reduction, and storage of calorimeter data are then

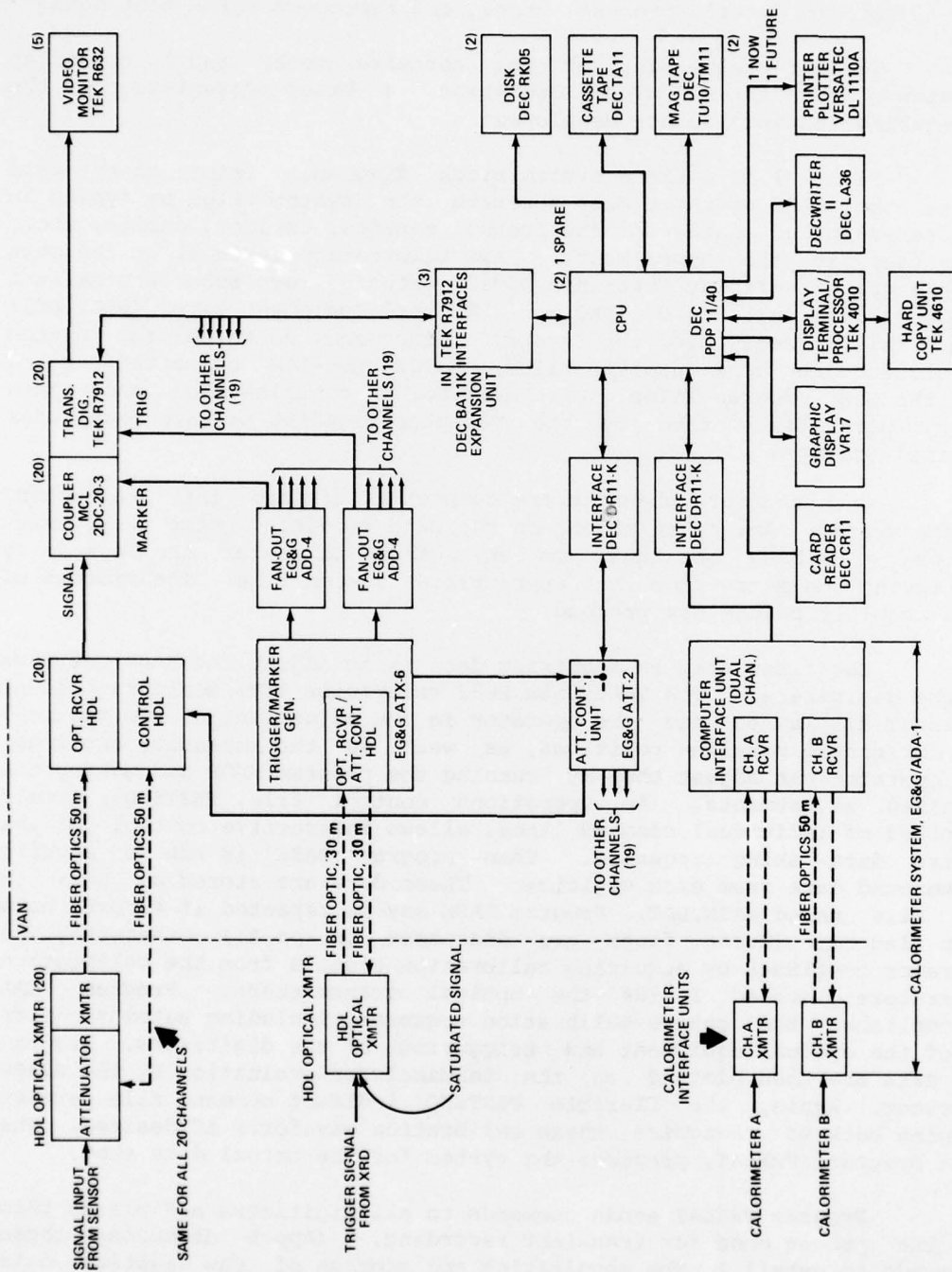


Figure 1. Systems generated electromagnetic pulse Transportable Automated Recording System (STARS).

performed by programs DATAPC, CONDAT, and MORDAT. Program DATAPC also turns on the optical transmitters and sets their programmable attenuators. As shown in figure 1, the interrupt trigger,  $\tau_0$ , is generated by an x-ray detector source and transmitted over a dedicated optical trigger link. It does so after program CONDAT has reached the wait mode signified by the message "WAITING FOR TRIGGER." The trigger, initiated by the simulator firing, is used to begin calorimeter data acquisition and to start the sweep for each digitizer. Sensor data are transmitted through the optical data links, where they are recorded by the digitizer's internal memory for a brief period before being transferred to DK1:.

Five steps remain:

- a. To normalize the data (sect. 4.1)
- b. To scale the data (program SCALE, sect. 3.2)
- c. To mark the data parameters (program MARK, sect. 4.2.4)
- d. To provide graphical and tabular output for the data user (sect. 4.3)
- e. To store the data on magnetic tape (sect. 4.4)

Section 5 includes a brief summary of all SASP programs. Detailed program descriptions may be found in the SASP Programmers' Handbook.<sup>2</sup>

## 2. DATA FILES

### 2.1 Introduction

System files are divided into three classifications: calibration, configuration, and data. Calibration files contain pertinent specifications on the measurement channel hardware--sensors, probes, optical data links, attenuators, baluns, etc. The system configuration file stores all the information needed to uniquely define the physical variables of the experiment at a given moment in time. In general, the SYSCON.DAT file is unique for any shot; however, if no variables are changed between shots and the simulator parameters remain unchanged,

---

<sup>2</sup>Robert Puttcamp, SKYNET Applications Software Package (SASP) Programmers' Handbook, Version 2, Harry Diamond Laboratories TR-1844 (February 1978).

only the shot number entry in the SYSCON.DAT file header will change. A copy of the current SYSCON.DAT file becomes an entry in each shot file. The shot file, Snnnnxx.DAT, in addition to the system configuration, contains all the data for a given shot. The data include the baseline, calibration, raw, and scaled data waveforms for each digitizer in the system, as well as computed parameters and user or operator comments. A complete description of all the files is presented in the following pages. A knowledge of these file structures is basic to understanding the operation of the software and usage of the data base.

Table I defines the data files used by SASP. The shot files and the MAIN.DAT file are on DK1.; the remaining are on DK0:. All are in user file directories (UFD's) with a user identification code (UIC) of [100,100].

TABLE I. DATA FILES

File	Data
PRBCAL.DAT	Probe calibration
ATNCAL.DAT	Device calibration
XCHCAL.DAT	Transmission channel calibration
SYSCON.DAT	System configuration
Snnnnxx.DAT (SHOT.DAT)	Taken for simulator shot named with four-digit shot number nnnn and two-character simulator code xx
MAIN.DAT	Preliminary name for Snnnnxx.DAT file

Once the shot number and simulator code locations in the SYSCON.DAT file have been set (by using program SHOTNO), program SYSCON points to the name of the current shot file. This name is composed of the letter "S" followed by the four-digit shot number and the two-character simulator code with the usual extension of "DAT". For example, for shot No. 115 on simulator P1, the shot file name is "S0115P1.DAT". This name is established for the file by program TABIN,

which in one of its functions renames the MAIN.DAT file to the shot file name pointed to by the SYSCON.DAT file. (The term "SHOT.DAT" is used in this handbook as a general reference to a shot data file.)

## 2.2 File Structures

In the following pages, the structure of the data files is laid out in both graphical and tabular form (tables I to IV, fig. 2 to 4). Wherever possible, space has been left for future expansion; these areas are designated as "spare."

Certain entries in the calibration files are for record keeping only and are not used by SASP. For the PRBCAL.DAT file, these entries are the following: measurement type, model, manufacturer's serial number, connector, and parameter; for the ATNCAL.DAT file: model and manufacturer's serial number; and for the XHCAL.DAT file: fiber number and receiver number. If delay calculations are not required, delay entries can be zero in the calibration files, as well as in the SYSCON.DAT file. The postreceiver gain value in each channel entry in the SYSCON.DAT file also is for record keeping. The manner of calculating the data link gain includes the data link gain in the system gain calculation automatically without needing an additional reference to this item.

The transmitters and channel numbers used in the SYSCON.DAT file are defined in the XHCAL.DAT file. Though the SYSCON.DAT file may have fewer channels defined than does the XHCAL.DAT file (never more), its first channel is the first in the XHCAL.DAT file and so on. To change transmitter addresses, the SYSCON.DAT file must be changed (by using program LOG). Changing digitizer addresses requires regenerating the SYSCON.DAT file by using program CONGEN.

The XRD bias item in the channel entries of the SYSCON.DAT file is a special record-keeping feature to be activated when the probe used is a biased x-ray detector (XRD).

The calibration files are handled by subroutines PRBREF, PRBSCH, ATNREF, ATNSCH, XCHREF, XCHDRI, and XCHSCH. The SYSCON.DAT file can be read and rewritten by the subroutine SCIO (SCIN, SCOUT as entry points). Subroutine FSCIN is used to read the SYSCON information from a shot file into common block SYSCON. The MAIN.DAT and shot files may be accessed from FORTRAN programs through the subroutine DATFIO. The calibration files are record oriented, while the SYSCON.DAT, MAIN.DAT, and shot files are block oriented.

TABLE II. CALIBRATION FILE ENTRIES

File	Size (bytes)	Item	Type <sup>1</sup>	Offset (bytes)
PRBCAL.DAT entry	4	Probe identification (four characters)	ASCII	0
	4	Transfer function (dB)	FP	4
	4	Delay	FP	8
	8	Units	ASCII	12
	2	Flag (such as integration required)	Integer	20
	2	Balun No.	Integer	22
	2	Measurement type	ASCII	24
	8	Model	ASCII	28
	8	Manufacturer's serial No.	ASCII	30
	8	Connector	ASCII	44
	4	Parameter (such as equivalent area)	FP	52
		54 Total		
		64 Provided for		
ATNCAL.DAT entry	4	Device identification	ASCII	0
	4	Gain (dB)	FP	4
	4	Delay	FP	8
	8	Model	ASCII	12
	8	Manufacturer's serial No.	ASCII	20
		28 Total		
		32 Provided for		
XCHCAL.DAT header	2	No. of channels	Integer	0
	4	$\Delta t_1$	FP	4
	4	$\Delta t_{trig1}$	FP	8
	4	$\Delta t_{trig2}$	FP	12
	4	$\Delta t_{fidu}$	FP	16
	4	$\Delta t_{vert}$	FP	20
	4	Balun No. 1 insertion gain	FP	24
	4	Balun No. 2 insertion gain	FP	28
	4	Balun No. 3 insertion gain	FP	32
		34 Total		
		48 Provided for		
XCHCAL.DAT entry	2	Transmitter No.	Integer	0
	2	Fiber No.	Integer	2
	2	Receiver No.	Integer	4
	4	Delay	FP	8
	4	Calibration pulse reference level	FP	12
		14 Total		
		16 Provided for		

<sup>1</sup>ASCII: American Standard Code for Information Interchange; FP: floating point.

TABLE III. SYSTEM CONFIGURATION (SYSCON.DAT) FILE ENTRIES

Entry	Size (bytes)	Item	Type <sup>1</sup>	Offset (bytes)
Header	6	Date	3 integers	0
	6	Time	3 integers	6
	4	Shot number	ASCII	12
	2	Shot number	Integer	16
	2	Number of digitizers	Integer	18
	2	Number of digitizer control boards	Integer	20
	2	Number of channels	Integer	22
	2	Virtual block number (VBN)	Integer	24
		baseline data		
	2	VBN calibration data	Integer	26
	2	VBN raw data	Integer	28
	2	VBN scaled data	Integer	30
	2	VBN parameters	Integer	32
	2	VBN user comments	Integer	34
	2	Simulator identification	ASCII	36
	2	Wire type code	ASCII	38
	4	Charging voltage output level	FP	40
	4	Tank vacuum	FP	44
	4	Voltage diode peak value	FP	48
	4	Current diode peak value	FP	52
	30	Comment (from TABIN)	ASCII	56
	2	Number of blocks in shot file	Integer	86
	4	Calorimeter 1 value	FP	88
	4	Calorimeter 2 value	FP	94
	2	Number of samples for calorimeter 1	Integer	96
	2	Number of samples, calorimeter 2	Integer	98
	2	Calorimeter error flag	Integer	100
	2	Filter code	Integer	102
	4	Calorimeter 1 uncertainty	FP	104
	4	Calorimeter 2 uncertainty	FP	108
	4	Filter thickness	FP	112
	2	Number of thicknesses of debris shield	Integer	116
	2	Debris shield code	Integer	118
	40	Spare		120
Digitizer block	2	R7912 device address	Integer	0
	2	Channel number	Integer	2
	2	Error flag	Integer	4
	2	Sweep rate	Integer	6
	2	Vertical sensor	Integer	8
	2	Bloom	Integer	10
	2	Baseline offset	Integer	12
	4	System gain in decibels	FP	14
	4	System gain factor	FP	18
	2	Spare		22
		(There are 20 of these blocks.)		
Channel block	2	Error flag	Integer	0
	6	Test point	RAD50	2
	4	Attenuator identification	ASCII	8
	4	Probe identification	ASCII	12
	2	XRD bias	Integer	16
	2	Transmitter control cable number	Integer	18
	2	Pretransmitter cable delay	FP	20
	4	Device identification	ASCII	24
	2	Transmitter number	Integer	28
	2	Transmitter attenuator setting	Integer	30
	4	Postreceiver cable delay	FP	32
	4	Postreceiver gain value	FP	36
	2	Number of digitizers on this channel	Integer	40
	2	Number of first digitizer	Integer	42
	2	Number of second digitizer	Integer	44
	2	Number of third digitizer	Integer	46
	4	Data link gain <sup>2</sup>	FP	48
	4	Data link scale factor <sup>2</sup>	FP	52
	4	Channel delay	FP	56
	4	Spare		60
		(There are 27 of these blocks.)		
Calibration Pointer	4	Calibration plus or minus value	FP	0
	4	Calibration zero value	FP	4
		(These two values are repeated for each digitizer in the system; 192 bytes are available for this information.)		

<sup>1</sup>ASCII: American Standard Code for Information Interchange.

FP: Floating point number.

RAD50: Base 50 representation of ASCII data.

<sup>2</sup>These values are not used; the real information is stored in the digitizer array. If there is more than one digitizer in the channel, the values for the last digitizer are stored here.

TABLE IV. MAIN DATA (SHOT.DAT) FILE ENTRIES

Entry	Size (blocks)	Item	Type <sup>1</sup>	Offset (bytes)
Components	5	Copy of finalized SYSCON		
	2 per TEK R7912	Baseline data		
	2 per TEK R7912	Calibration data		
	8 per TEK R7912	Raw data		
	4 per TEK R7912	Scaled data		
	1/4 per TEK R7912	Parameters		
	1/8 per TEK R7912	User comments (64 characters) last block filled		
	(bytes)			
Parameters	2	Evaluation flag	Integer	0
	2	Fiducial time	Integer	2
	2	Start of data time	Integer	4
	2	Time of peak of interest	Integer	6
	2	End of data time	Integer	8
	2	Time of peak of calibration	Integer	10
	4	Time to start of data from fidu	FP	12
	4	Peak of interest	FP	16
	4	Time to peak of interest from data start	FP	20
	4	Area	FP	24
	24	Comment on calibration data	ASCII	28
	24	Comment on raw data	ASCII	52
	24	Comment on scaled data	ASCII	76
	100	Total		
	128	Provided for		

<sup>1</sup>FP: Floating point;

ASCII: American Standard Code for Information Interchange

PRBCAL.DAT

RECORD NO. 1	FIRST PROBE ENTRY	32-WORD RECORD
	⋮	
	LAST PROBE ENTRY	32-WORD RECORD
	-1'S	32-WORD RECORD

ATNCAL.DAT

RECORD NO. 1	FIRST ENTRY	16-WORD RECORD
	⋮	
	LAST ENTRY	16-WORD RECORD
	-1'S	16-WORD RECORD


XCHCAL.DAT

RECORD NO. 1	HEADER, PART I	8-WORD RECORD
RECORD NO. 2	HEADER, PART II	8-WORD RECORD
RECORD NO. 3	HEADER, PART III	8-WORD RECORD
RECORD NO. 4	FIRST CHANNEL ENTRY	8-WORD RECORD
	⋮	
RECORD NO. (3+m)	LAST CHANNEL ENTRY	8-WORD RECORD

NOTE: m IS NUMBER OF CHANNELS AS DEFINED IN HEADER.

Figure 2. Calibration file structures.

Figure 2. Calibration file structures (cont'd).

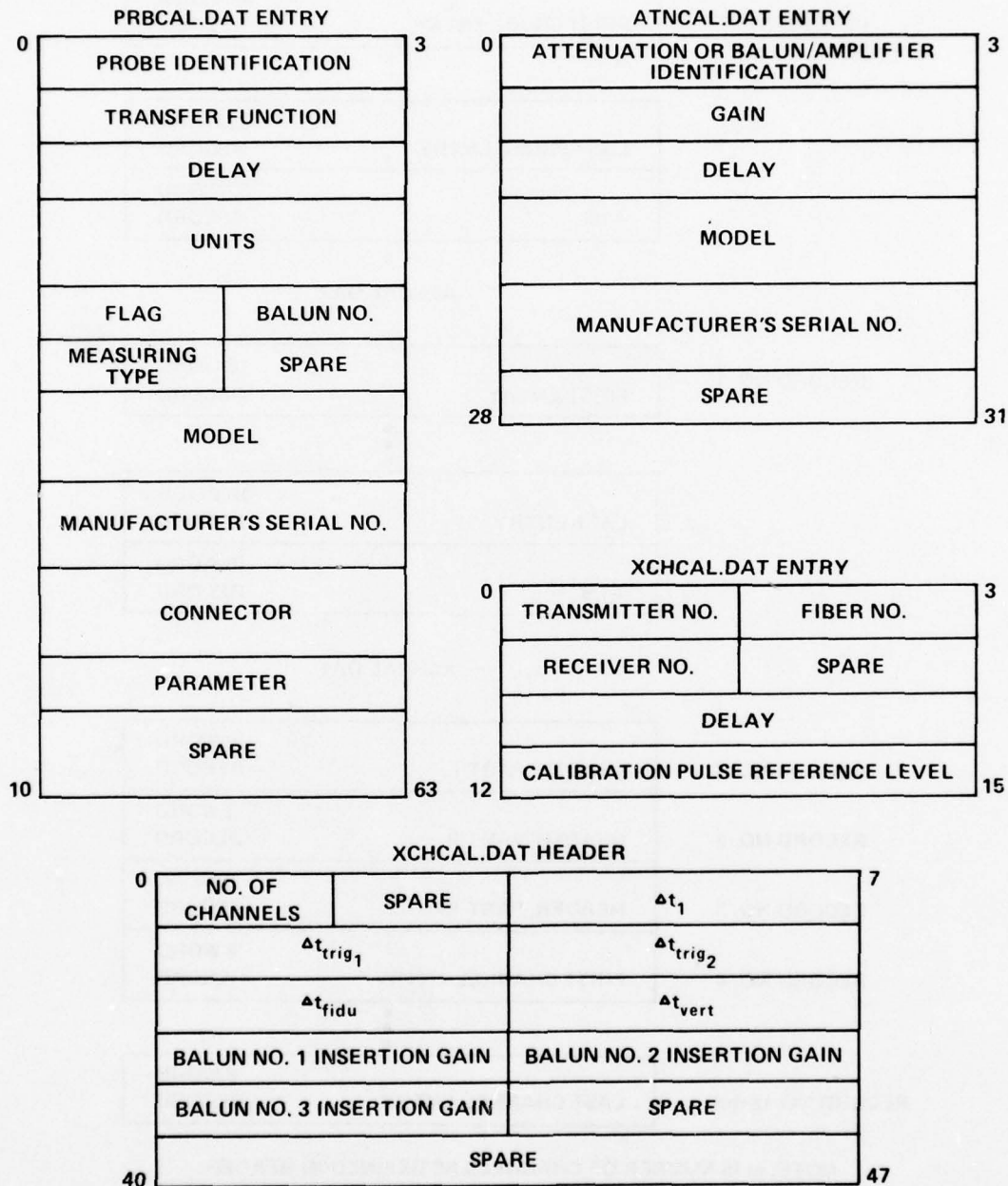

NOTE: NUMBERS OUTSIDE OF BOXES ARE BYTE NUMBERS.


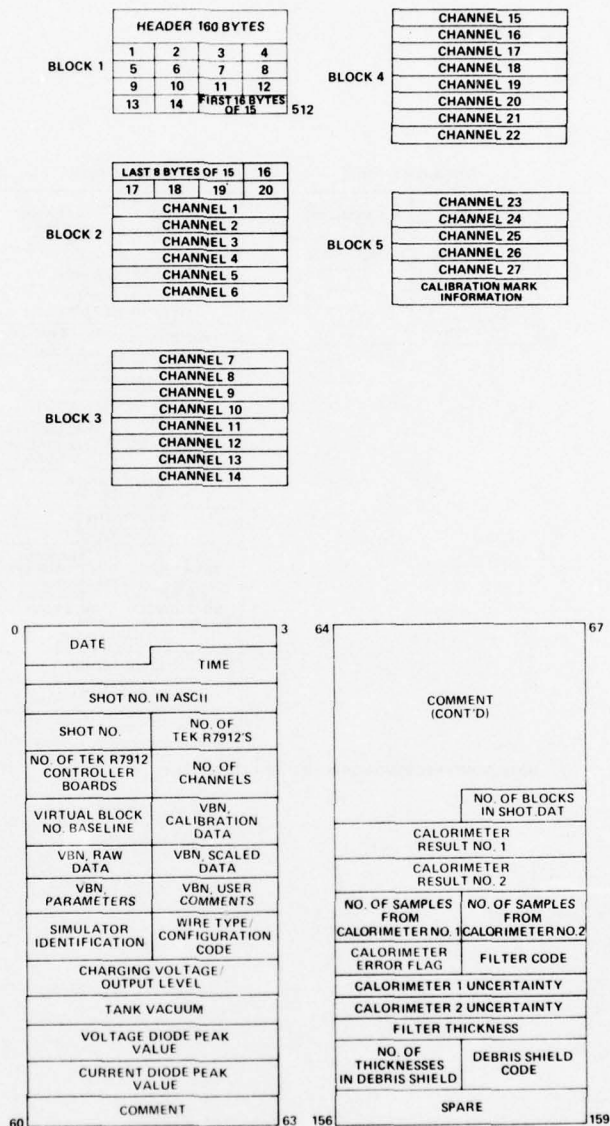


Figure 3. System configuration (SYSCON.DAT) file structure.

Figure 3. System configuration (SYSCON.DAT) file structure (cont'd).

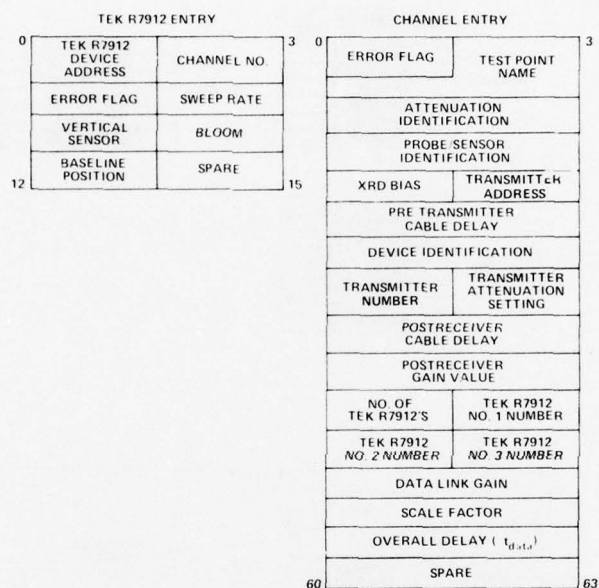

NOTE: NUMBERS OUTSIDE OF BOXES ARE BYTE NUMBERS.

# COMPONENTS


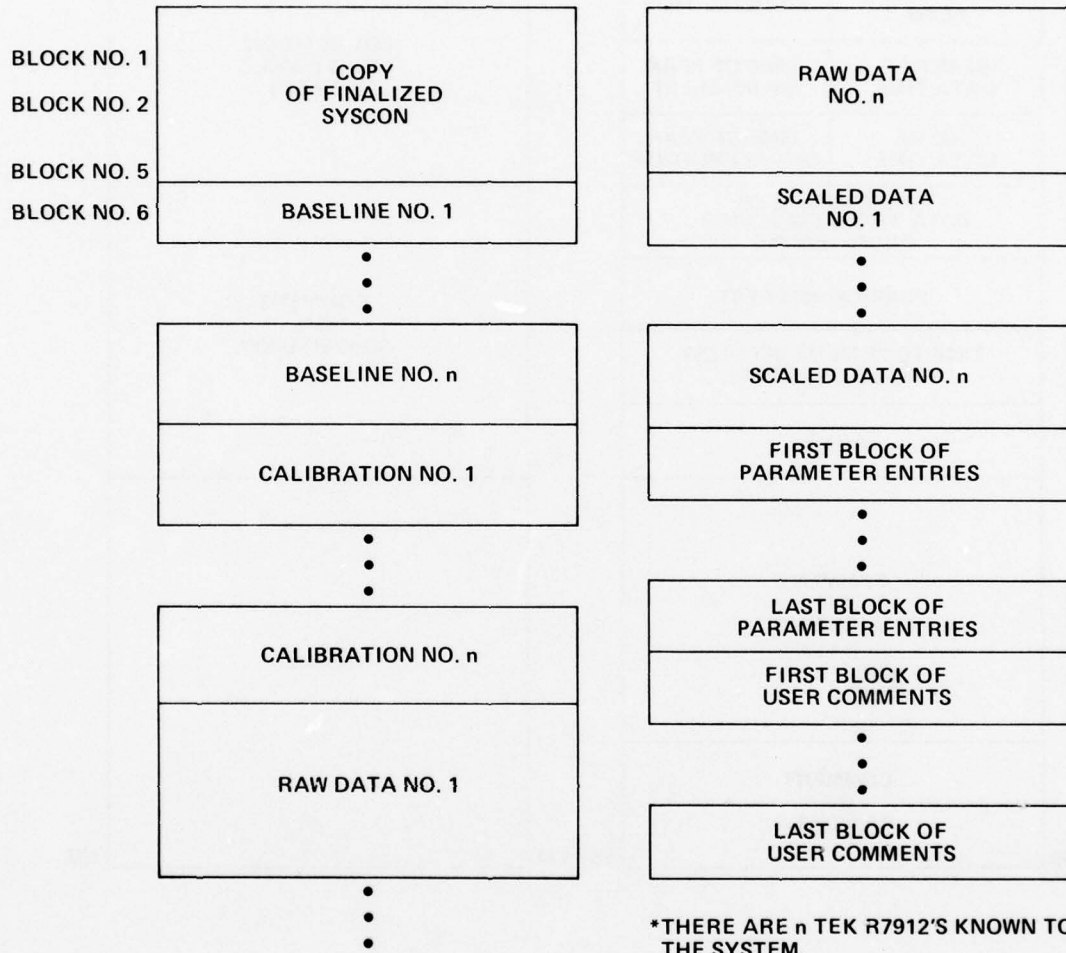


Figure 4. Shot file structure.

Figure 4. Shot file structure (cont'd).

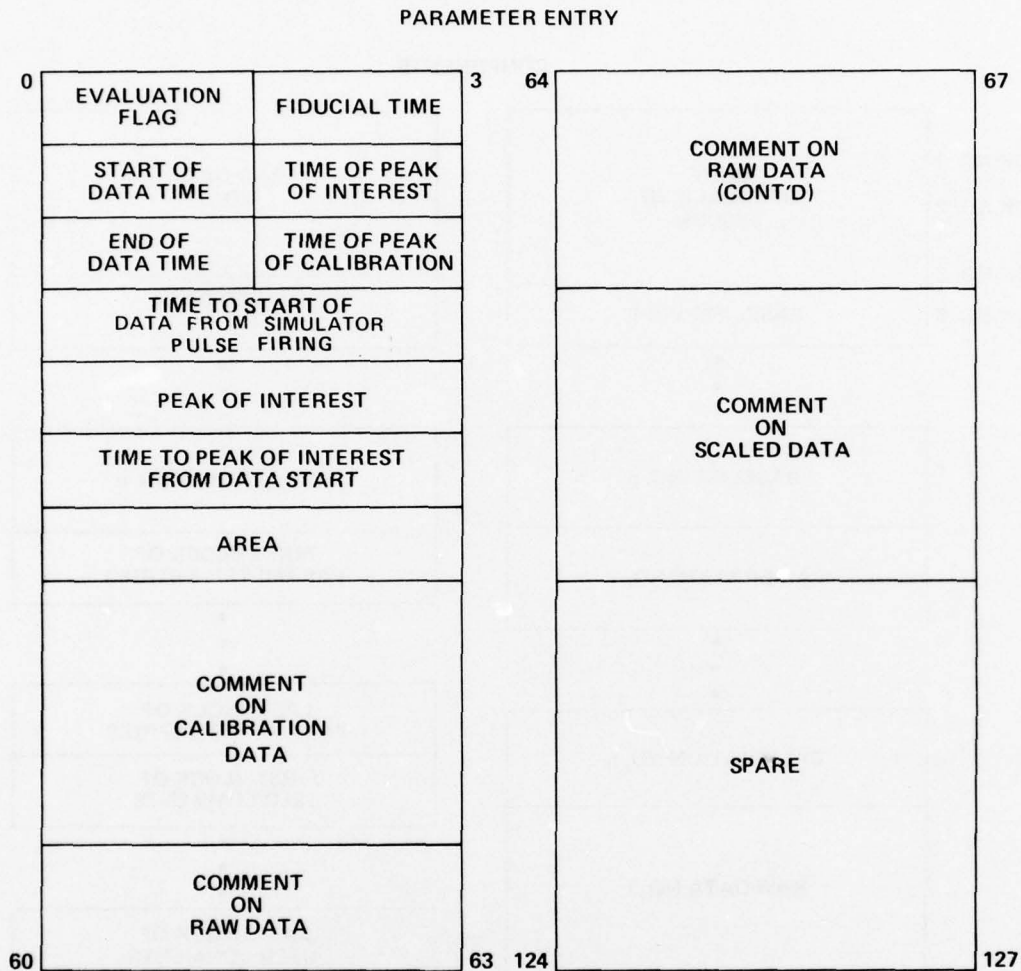

### 3. DATA SYSTEM OPERATIONS

#### 3.1 Introduction

Prior to testing, many unknown quantities remain to be defined. The data system calibration files must be generated; information concerning physical test variables must be stored in memory. Amplifier

gains, calibration pulse reference levels, and sensor and probe identification numbers must be made known to the data system. Section 3.2 defines the procedures for constructing calibration files, the files that contain the information necessary to conduct a test. Separate programs exist for modifying each of these calibration files. Log forms to be filled out by the STARS operator and the experimenter also are described. The information on these logs provides the needed cross references between calibration and data files.

Next, terminal input conventions are presented. This information is of use to the STARS console operator for the input of tabular data from the log forms.

A description of the data-taking sequence follows in section 3.4. The SASP software can be run in the manual or automatic mode. That is, the separate, modular programs can be combined into an indirect command file and run automatically. The SASP software was specifically designed to run in this mode, under the control of TESTSEQ. Section 3.4.1 explains the features of TESTSEQ and the procedures for abnormally aborting TESTSEQ. Precautions to be observed for running the SASP programs in the manual mode also are given. Discussion on disc-to-tape transfer conventions and the method for preparing new discs conclude this section.

### 3.2 Data System Generation

The calibration and configuration files, XCHCAL.DAT, ATNCAL.DAT, PRBCAL.DAT, and SYSCON.DAT, must be established before data taking can occur. Devices and attenuators to be used should be entered in the ATNCAL.DAT file by using program ATN. Probes are specified in the PRBCAL.DAT file with program PRB.

The contents of the XCHCAL.DAT file must be established before generating the SYSCON.DAT file; otherwise, the transmitter numbers are not available when needed. Program XCH provides for the generation of the XCHCAL.DAT file. Then program CONGEN may be used to establish the SYSCON.DAT file.

Before beginning data system generation, the alignment of transmitters, channels, and digitizers should be decided. Creating the XCHCAL.DAT file provides the association between the channel number and the transmitter numbers. After the SYSCON.DAT file is created, program LOG is run, and the selected channel numbers are input into the digitizer log. These actions complete the link.

To change a digitizer address, the SYSCON.DAT file is regenerated. To modify a transmitter address, program LOG must be run. The SASP data input requirements follow.

The generalized SASP data channel can be visualized as shown in figure 5. Each item has a gain or delay or both associated with it, Gains are given in decibels calculated from equations of this form:

$$\text{Attenuator gain (dB)} = 20 \log_{10} \frac{\text{output (selected units)}}{\text{input (same units)}} ,$$

$$\text{Probe gain (dB)} = 20 \log_{10} \frac{\text{output (volts)}}{\text{input (selected units)}} ,$$

$$\text{Device gain (dB)} = 20 \log_{10} \frac{\text{output (volts)}}{\text{input (volts)}} ,$$

The gain of the transmission channel ("data link gain") is calculated from a knowledge of the calibration pulse reference level and a measurement of the peak recorded at the digitizer. The gain is calculated in program CHKCAL (sect. 4.2.2). All the gains are combined into a single overall scale factor by program SCALE. There, the overall delay also is calculated.

The Tektronix R7912's have many characteristics. They receive their own entry in the configuration file (SYSCON.DAT) which is linked to the channel entry by the channel number. Three digitizers may be so linked to a single channel. (See the SASP file structures, sect. 2.2.)

The attenuator, probe, device, and transmission channel also receive special treatment. Calibration files are established to define their characteristic gains and delays and to perform other bookkeeping functions. The input required is as follows:

Device or attenuator (ATNCAL.DAT file, program ATN)

Identification (four characters)  
Gain (decibels)  
Delay (nanoseconds)  
Model (eight characters)  
Manufacturer's serial number (eight characters)

Probe (PRBCAL.DAT file, task PRB)

Identification (four characters)  
Transfer function (gain in decibels)  
Delay (nanoseconds)  
Data processing flag  
1 = integration  
0 = simple scaling

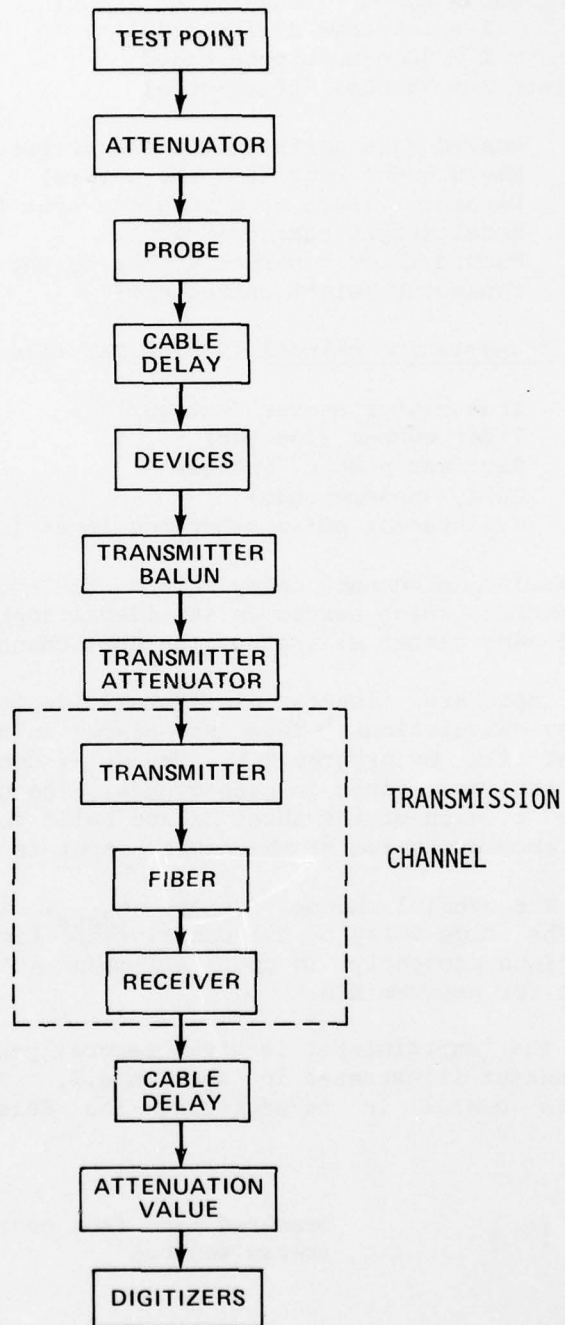


Figure 5. Generalized data channel.

Balun number (transmitter input)

1 = 100 ohms differential

2 = 50 ohms single ended

3 = 78 ohms differential

Scaled data units (eight characters) (selected units)

Measurement type (two characters)

Parameter (such as equivalent area for B-dot probe)

Model (eight characters)

Manufacturer's serial number (eight characters)

Connector (eight characters)

Transmission channel (XHCAL.DAT file, program XCH)

Transmitter number (integer)

Fiber number (integer)

Receiver number (integer)

Delay (nanoseconds)

Calibration pulse reference level (millivolts)

Each transmission channel entry in the XHCAL.DAT file is assigned a channel number, which serves as its identification and which is forced to be the same number as that of the data channel that includes it.

There are also certain system-wide constants required for gain and delay calculations. These are placed in the header section of the XHCAL.DAT file by program XCH. The delay constants are illustrated in figure 6 and are given in nanoseconds. The gain values are the insertion loss of each of the three baluns built into the transmitters. The negative absolute value of the number input is stored.

The overall channel delay,  $\Delta t_{data}$ , is calculated in program SCALE. The time delay of the oscilloscope time base,  $\Delta t_e$ , drops out of any equations presently in use. All other  $\Delta t$ 's illustrated in figure 6 are input for program XCH.

The experimenter is given several pertinent delay times on the tabular output illustrated in section 4.3. The following definitions will prove useful in interpreting the delay time factors in these tables:

DT1 ( $\Delta t_1$ )

measured time from source to trigger sensor  
(x-ray sensor)

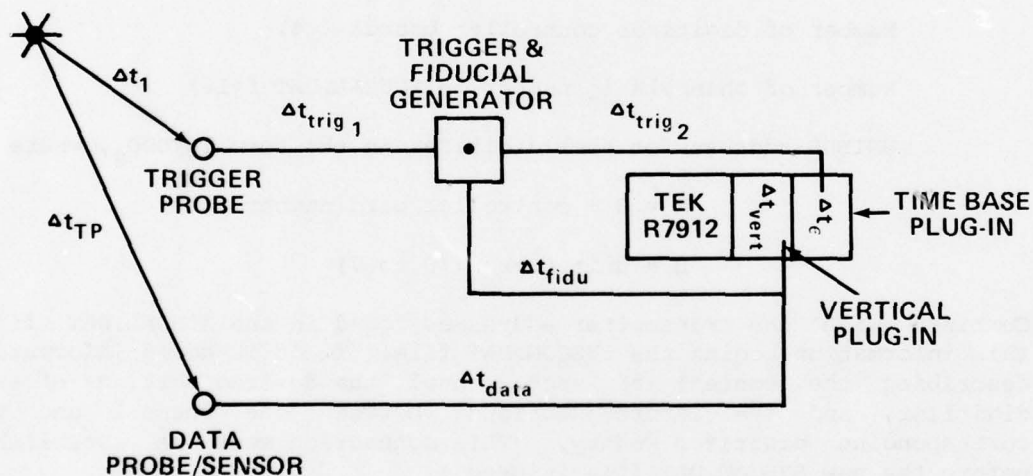
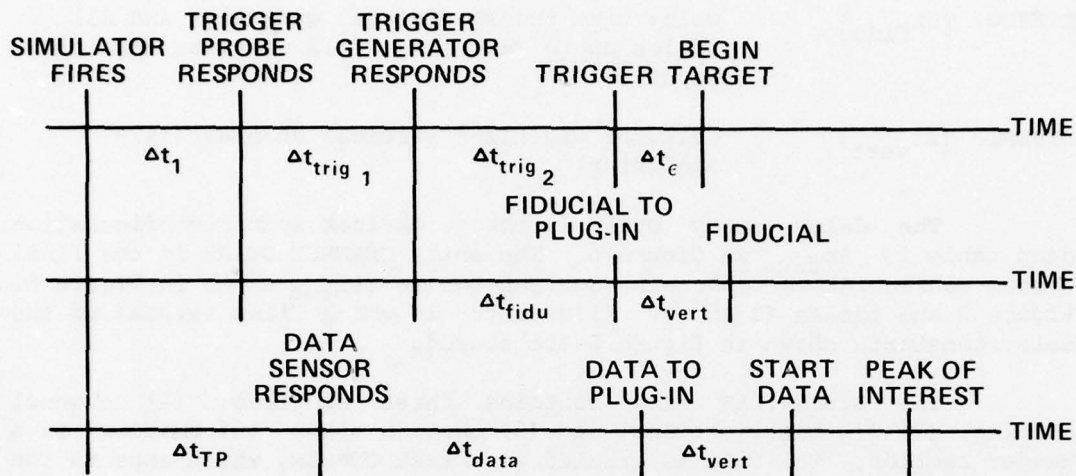


Figure 6. Delay considerations.

DTRIG2 ( $\Delta t_{\text{trig}_1} + \Delta t_{\text{trig}_2}$ )

trigger sensor delay from trigger sensor to fiber optic transmitter, fiber optic trigger link delay, delay time through trigger generator, and delay in all cables from output of trigger generator to external trigger input on Tektronix R7912

DTFIDU	$(\Delta t_{\text{fidu}})$	delay time through fiducial generator and all cables up to Tektronix R7912 vertical amplifier input
DTVERT	$(\Delta t_{\text{vert}})$	delay of digitizer vertical channel (7A19 amplifier)

The delay entry DELAY in the quick-look system configuration dump table is  $\Delta t_{\text{data}}$  in figure 6. The entry CHANNEL DELAY in the final system configuration table also corresponds to  $\Delta t_{\text{data}}$  shown in figure 6. Figure 2 and tables II to IV illustrate in which files several of the delay constants shown in figure 6 are stored.

The SYSCON.DAT file contains three sections: (1) channel entries, (2) digitizer entries, and (3) certain other information in a header section. The file is created with task CONGEN, which expects the following input:

Number of digitizers ( $\leq 20$ )

Number of digitizer controller boards ( $\leq 4$ )

Number of channels ( $\leq$  number in XHCAL.DAT file)

UNIBUS address for each digitizer in the form  $16X0D0_8$ , where

$X = 3 + \text{controller card number}$

$D = \text{unit number (0 to 7)}$

Combined with the transmitter addresses found in the XHCAL.DAT file, this information begins the SYSCON.DAT file. To it is added information describing the content of each channel, the desired settings of each digitizer, and the interconnections between the channel and the corresponding digitizer entry. (This connection must be established before the new SYSCON.DAT file is used.)

Four input log forms are necessary to fulfill the entry requirements of the SYSCON.DAT file. The information is entered with program LOG, which has four "pages" of input that correspond to the format of the logs. The standard SASP terminal input conventions apply.

Figures 7 to 10 are sample logs. The comments column on the transmitter log (fig. 8) is for user comments to be entered at a different time into a different file. The exact nature of the entries is as follows:

Probe log (fig. 7)

TEST PT: test point name, nine alphanumeric characters maximum for specifying probe location, test conditions, and other characteristics

ATTN ID: identification of attenuator entered in ATNCAL.DAT file ("blank" means no attenuator)

SENSOR ID: identification of probe entered in PRBCAL.DAT file

PRE XMTR CABLE DELAY: value to reflect delays not accounted for by other pretransmitter entries

DEVICE ID: identification of device entered in ATNCAL.DAT file ("blank" means no device)

PROBE LOG

SHOT NO. \_\_\_\_\_ (next expected)

CHANNEL	TEST PT	ATTN ID	SENSOR ID	PRE XMTR CABLE DELAY	DEVICE ID
1					
2					
3					
.					
.					
.					
20					

Figure 7. Probe log.

Transmitter log (fig. 8)

TRANSMITTER ATTN: transmitter attenuation, 0 to 45 dB in 3-dB steps (negative absolute value of input is stored)

XRD BIAS: special case, integer parameter to be used in specifying variable probes, particularly biased XRD's

TRANSMITTER ADDRESS: control cable number from EG&G attenuator controller box

COMMENTS (LAST SHOT): user comments

#### TRANSMITTER LOG

SHOT NO. \_\_\_\_\_ (next expected)

CHANNEL	XMTR ATTN	XRD BIAS	XMTR ADDRESS	COMMENTS (LAST SHOT)
1				
2				
3				
.				
.				
.				
20				

Figure 8. Transmitter log.

#### Receiver log (fig. 9)

DELAY (REC-DIG): value reflecting delay between receiver and digitizer

ATTN (REC-DIG): value recording attenuation between receiver and digitizer (not used in calculations, since measuring data link gain automatically includes it)

# RECEIVER LOG

DATE \_\_\_\_\_

START SHOT \_\_\_\_\_ END SHOT \_\_\_\_\_

CHANNEL	DELAY (REC-DIG)	ATTN (REC-DIG)
1		
2		
3		
.		
.		
.		
20		

Figure 9. Receiver log.

## Digitizer log (fig. 10)

DIG. NO. (OPER): operator's alignment of channels and digitizers; program LOG prints digitizer number and requests channel number

SWEEP (NS/DIV): sweep rate in nanoseconds/division

VERT SENS (MV/DIV): vertical sensitivity in millivolts/division

BASELINE OFFSET ( $\pm$  MINOR DIV): must be zero---this feature is not currently implemented

CHANNEL NO.: channel to which this digitizer is connected

Simulator information is entered into the header area of the SYSCON.DAT file in two stages. Program SHOTNO requests and stores a one-to four-digit shot number and a two-character simulator code, which should be "O2" for the Oil Water Line II, SGEMP Simulator (OWL II). The

# DIGITIZER LOG

DIG. NO. (OPER)	SWEEP (NS/DIV)	VERT SENS (mV/DIV)	BASELINE OFFSET ( $\pm$ MINOR DIV)	CHANNEL NO.

Figure 10. Digitizer log.

second stage takes place with program TABIN, which has two modes of input, a special one (for O2) and a general one.

## Special (O2) mode

Wire type (two characters)  
Pulse charge (megavolt)  
Tank vacuum (torr)  
Voltage diode peak (volt)  
Current diode peak (volt)  
Comment (30 characters)  
Debris shield code (integer)  
Number of debris shield thicknesses (integer)  
Filter code (integer)  
Filter thickness (floating point)

## General mode

Configuration code (two characters)  
Output level (units not yet established)  
Comment (30 characters)

Space in the shot file has been provided for user comments. They are logged on the transmitter log for the shot to follow. The program used to enter them is USRCOM. It expects 64 or fewer characters.

In all terminal entries, except the digitizer UNIBUS address in input to the CONGEN file, standard SASP terminal input conventions apply.

### 3.3 Terminal Input Conventions

A line of input is read in by subroutine INPUT. This line consists of fields (or columns) separated by the control I or TAB character. The first field begins with the first character. Each subsequent control I begins a new field, irrespective of preprinted column headings.

The three forms of input conversion--American Standard Code for Information Interchange (ASCII), integer, and floating point--are performed by subroutines ASDCD, FXDCD, and FLDCD, respectively. The choice of field is made by a field-number parameter. See the individual documentation for each subroutine.<sup>2</sup> The DEC terminal input line editing commands, control U and RUBOUT, are still operative.

For a given field, there are three null conditions:

- a. The input was terminated (carriage return) before the field was reached.
- b. The field was reached, but a second control I was given immediately.
- c. The first character of the field was a blank.

Under conditions a and b, the value stored before input remains untouched. Condition c causes a zero to be stored.

### 3.4 Data-Taking Sequence

#### 3.4.1 TESTSEQ Operations Control File

The SASP data-taking sequence is controlled by indirect command file DK1:[100,100]TESTSEQ.CMD. (For general information on indirect command files, see the Operator's Procedures Manual for DEC RSX-11M.<sup>3</sup>) Decisions by the operator are entered by means of yes (Y) or no (N) answers to preprogrammed questions. Operations are carried out by programs executing on command from TESTSEQ. Many programs require operator input or provide output. Pauses in command file processing are placed at crucial points; the system is then under the direct control of the operator.

---

<sup>2</sup>Robert Puttcamp, SKYNET Applications Software Package (SASP) Programmers' Handbook, Version 2, Harry Diamond Laboratories TR-1844 (February 1978).

<sup>3</sup>Operator's Procedures Manual RSX-11M Version 2, DEC-11-OMOGA-B-D, Digital Equipment Corp., Maynard, MA (September 1975).

If abnormal conditions require exiting from TESTSEQ's control, typing control C and a carriage return "punches through" to the operating system. If TESTSEQ has paused, is waiting for an answer to a question, or is printing comments, the operator types "ABO AT. (CR)". If a SASP program is running, the operator should type "ABO AT. (CR)" and "ABO" plus first three characters of the program name. This will abort both the running program and the indirect command file processor. If only the command "ABO AT." is entered, the indirect command file processor will be aborted at the termination of the running program. If the aborted program referenced a data file, the file may be locked. As soon as possible, the operator should use the DEC utility PIP (Peripheral Interchange Program) to unlock the file.

Any task that references the digitizers or the transmitters may produce error messages at the terminal. The standard SASP digitizer and channel error codes apply (see app A for details).

Before the data-taking sequence can begin, a few operations must be performed or confirmed. The RSX-11M operating system must be running from a suitable system disc: one that has all the needed DEC utility tasks and SASP data files. A data disc with sufficient storage space must be mounted; if it is not, the operator makes room or initializes a new one. The correct time and date should be entered with the monitor console routine (MCR) command TIME. The indirect command file DATAOP.DAT should be processed, either as a direct MCR command or by use of TESTSEQ.

Shot files are to be transferred to magnetic tape as quickly as possible. Also, at some established interval (such as once a day), the master data tapes should be updated.

Figure 11 is the flowchart for TESTSEQ. Errors during task BASE should be resolved before TESTSEQ continues. Besides acquiring baselines, this task is meant to act as a screen for all potential problems at other, more crucial steps in the process. The TESTSEQ Operations Control File (app C) can be used as an aid to understanding the flowchart. Individual task documentation in the SASP Programmer's Handbook,<sup>2</sup> as well as the source listings, also will prove helpful.

---

<sup>2</sup>Robert Puttcamp, *SKYNET Applications Software Package (SASP) Programmers' Handbook, Version 2*, Harry Diamond Laboratories TR-1844 (February 1978).

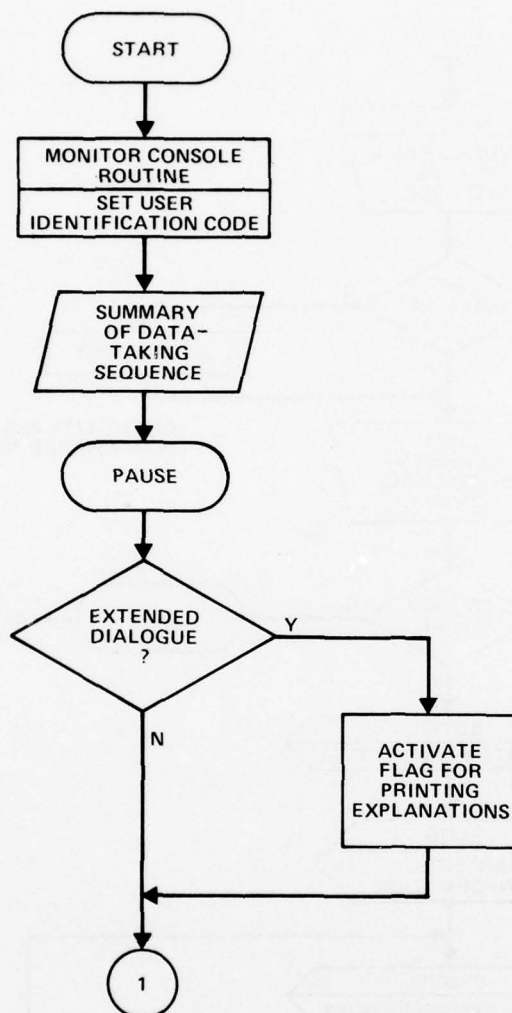


Figure 11. Flowchart for TESTSEQ.

Figure 11. Flowchart for TESTSEQ (cont'd).

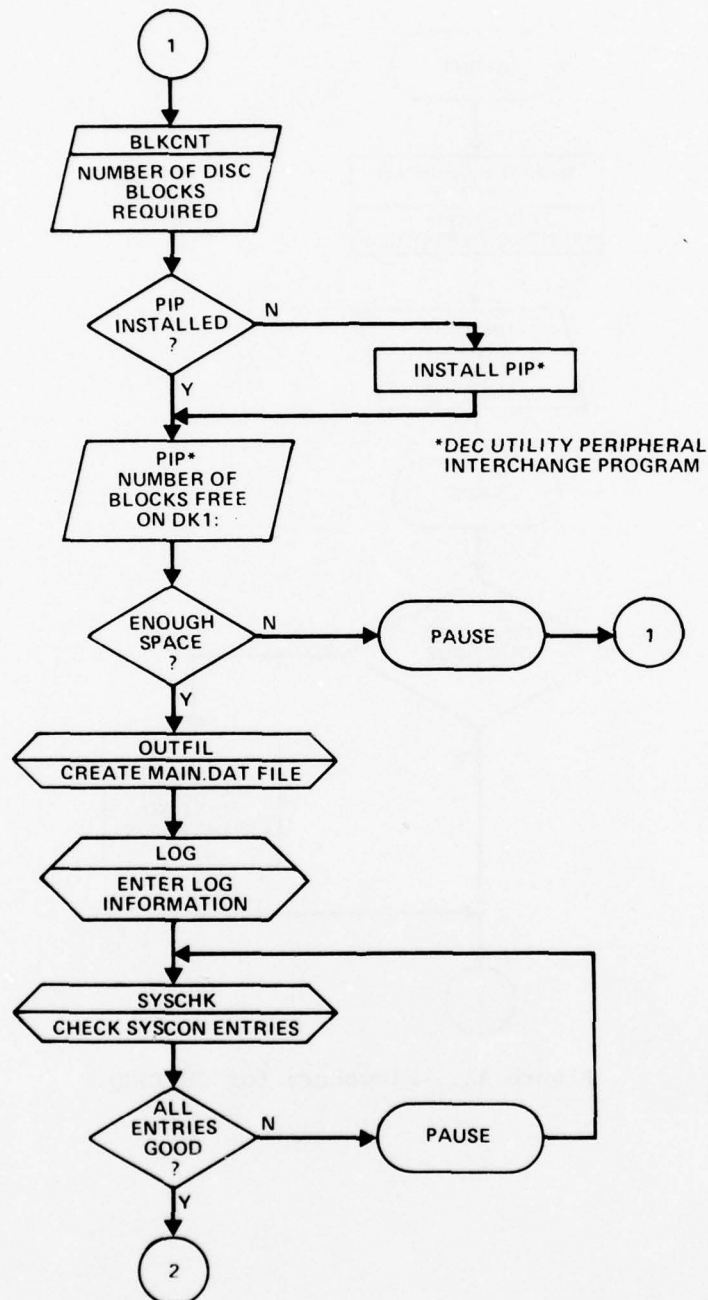


Figure 11. Flowchart for TESTSEQ (cont'd).

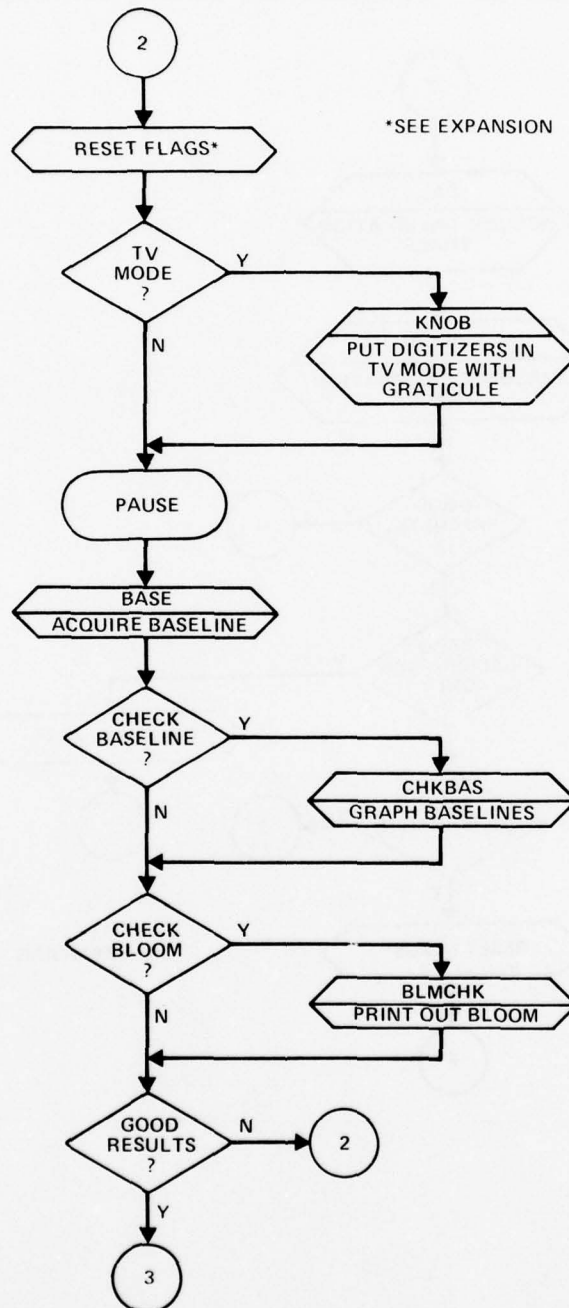


Figure 11. Flowchart for TESTSEQ (cont'd).

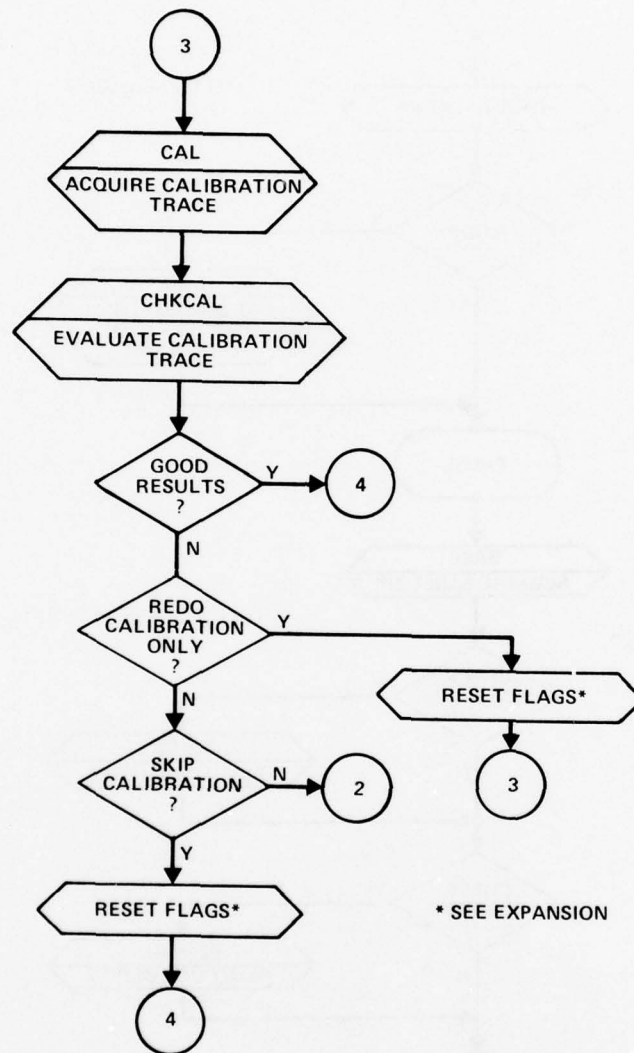


Figure 11. Flowchart for TESTSEQ (cont'd).

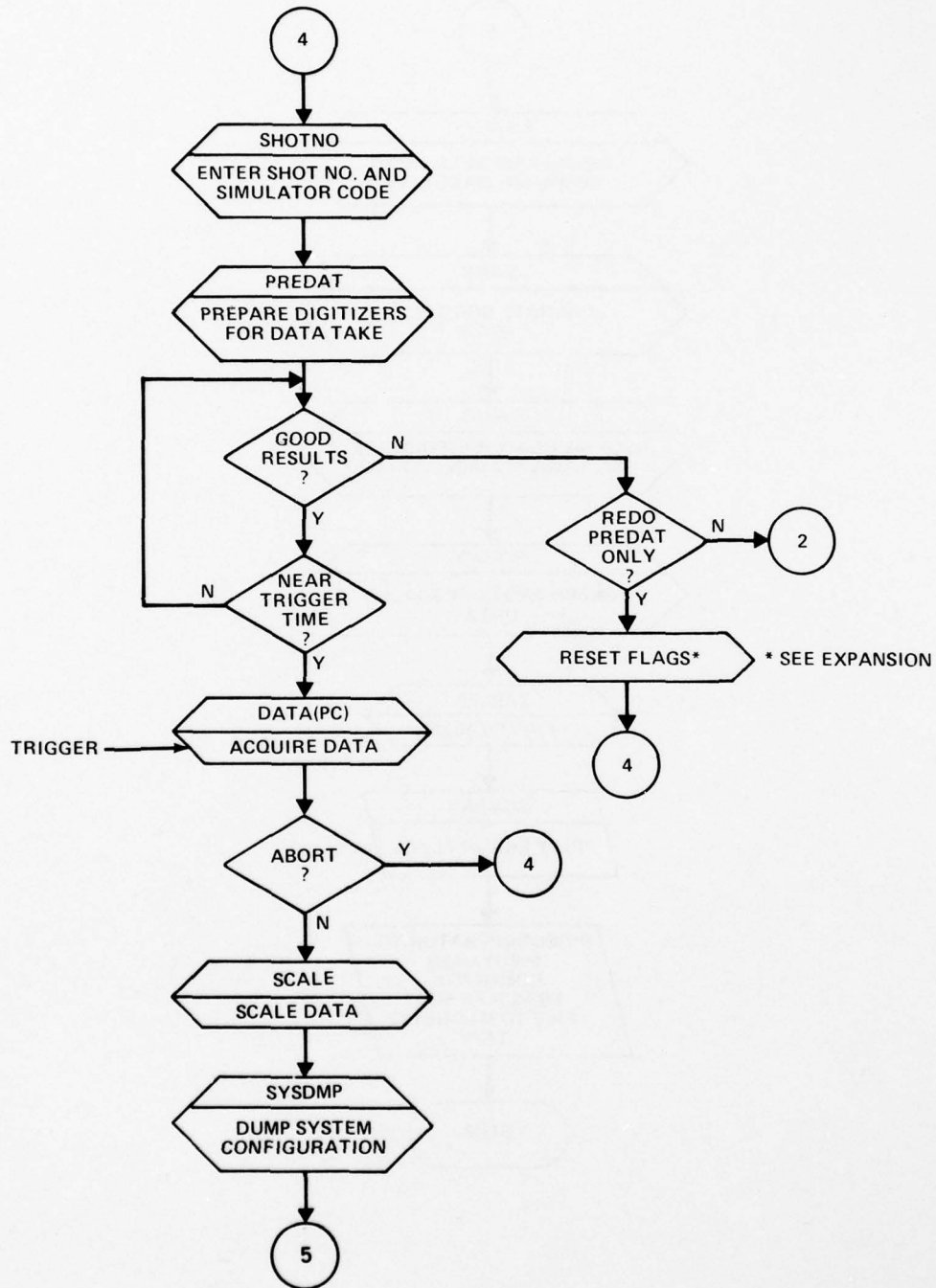


Figure 11. Flowchart for TESTSEQ (cont'd).


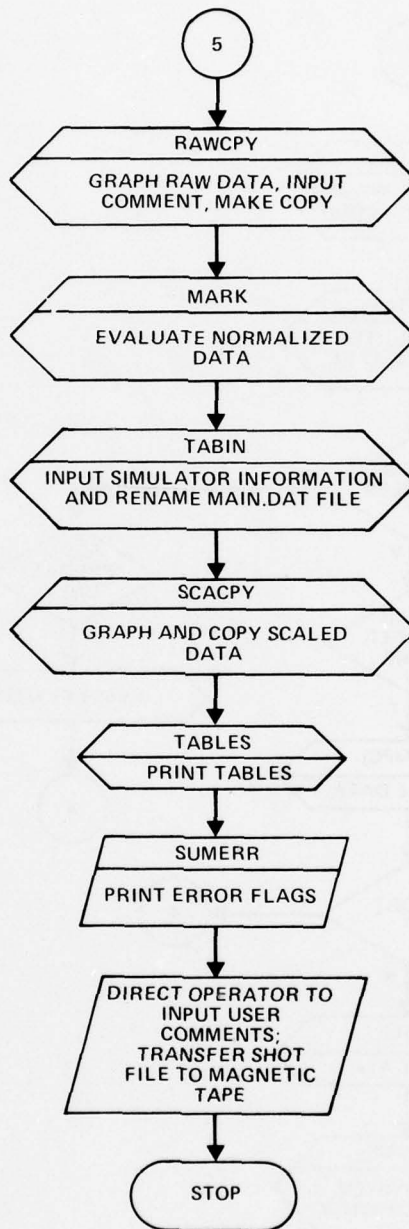


Figure 11. Flowchart for TESTSEQ (cont'd).

EXPANSION OF

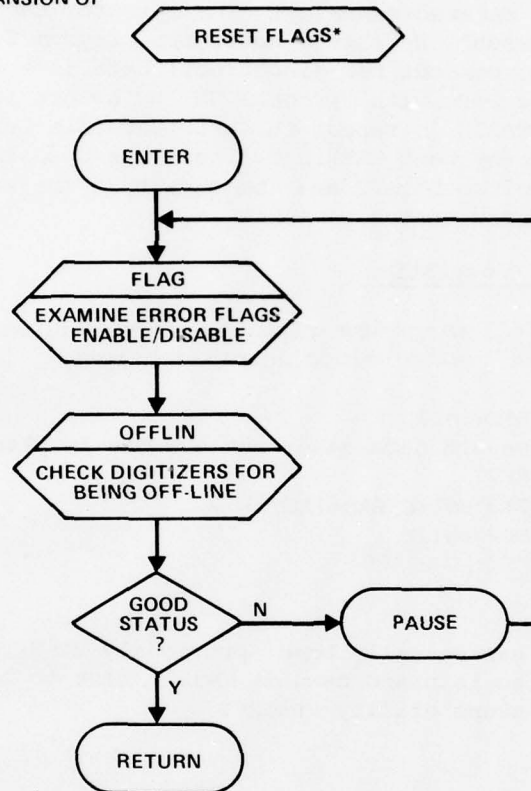

### 3.4.2 Manual Data Acquisition Procedures

Care must be taken when TESTSEQ is abandoned and data-taking instructions are entered manually. The operator should be very familiar with the processes. He should never run program SCALE without running task CHKCAL to screen out bad (incorrect) data link gains. The operator sets the shot number with program SHOTNO before producing any output. He runs program TABIN to rename the MAIN.DAT file before running program OUTFIL to create the next MAIN.DAT file. The operator should always run program PREDAT before DATAPC and be sure that the error flags represent a desired and proper status.

### 3.4.3 Disc Preparation

The following summarizes the instructions for preparing a new data disc ("valid" means volume identification):

```
DMO DK2:valid
(Remove old data disc; put new one in place.)
BAD DK2:
INI DK2:valid/BAD=[AUTO]
MOU DK2:valid
UFD DK2:[100,100]
```

If the new disc has actually been previously used, the contents will be lost. If the disc is brand new, it will have to be formatted first by using the stand-alone utility PRESRV.

## 4. DATA PROCESSING

### 4.1 Normalization

The raw data acquired through digitizing a waveform with the digitizer are not in a form compatible for processing by a digital computer. Usually, there is an overabundance of data describing the waveform. Sometimes, there is a lack of data describing certain portions of the waveform. Ideally, there are two vertical values for each horizontal value. These vertical values correspond to the upper and lower edges of the path traced by an electron beam as it sweeps across the target. Here, a simple average can determine the true center of the trace to a resolution of 1 part in 512. For more or less than two vertical values, an algorithm is required to determine the most likely center of the trace for that particular point on the raw data waveform. It is the task of the normalization routines to process the

raw data, with their excesses and deficiencies, and to obtain from them a more usable presentation of the waveform.

There are six basic steps to the normalization:

- a. Average to the center of the trace (ATC).
- b. Eliminate noise points.
- c. Fill gaps by linear interpolation.
- d. Fill end gaps with zeros.
- e. Subtract the baseline trace.
- f. Subtract 256.

For baseline traces and calibration traces, these steps are carried out with integer arithmetic. Data-trace normalization uses floating point operations. Step e is not used for baseline traces; in step d, a baseline trace is extrapolated rather than filled with zeros.

a. The ATC reduces the raw digitizer data to, at most, a single vertical value for each horizontal value. The original raw data are composed of horizontal and vertical values in a one-dimensional array. The vertical values that are associated with a horizontal value follow it. Each such group is in order by horizontal value, except that the data start at an arbitrary horizontal position of the target and wrap to the starting point. Horizontal values are differentiated from vertical ones by having the bit in variable HZMSK set (app D). ("HZMSK" is defined as  $1000_8$  by the SYMDEF macro, a macro used to initialize all SASP parameters. See app D.) Points flagged with the bit in variable BADBIT are excluded from consideration. ("BADBIT" is defined as  $10000_8$  by the SYMDEF macro.) If a horizontal value is flagged as bad, it is assumed that the associated vertical values also are bad. Vertical values are recorded from the top of the target to the bottom.

Once a vertical value is masked with variable DTAMSK, its value is the vertical position on the target and is in the range of 0 to 511. When a horizontal value is masked, the result is the complement of the position on the target. The complement also is in the range of 0 to 511. The coordinates (0,0) correspond to the lower left-hand corner. Each horizontal-vertical group is resolved into a horizontal value with a single vertical value.

The ATC routines consider four groups: (1) a single vertical value, (2) two vertical values, (3) three vertical values, and (4) four or more vertical values. Each group is considered independently.

Group 1.--It is assumed that the bottom edge of the trace was missed due to the dead time of the read gun; the dead time is equivalent to the time for the scan gun to read four points. The single vertical value is then shifted down slightly by subtracting two.

Group 2.--The separation of the two vertical values is compared to variable BLMFAC times the bloom parameter. If the separation is less than this product, the two are averaged. If not, the separation is flagged as a gap point.

Group 3.--The closest pair is chosen, and the conditions of group 2 are applied.

Group 4.--Four or more vertical values are immediately flagged as a gap point. The flag for a gap point is "-1." All well-defined vertical values are from 0 to 511.

Since the bloom parameter used in groups 2 and 3 is calculated from the baseline trace (see documentation for subroutine BASATC in the Programmers' Handbook<sup>2</sup>), the baseline ATC requires a simplified approach. Group 1 remains the same. For group 2, the averaging takes place automatically. Groups 3 and 4 result in gap points.

b. The trace is now a single dimensional array of 512 vertical values and gap points with the horizontal value being implied by the subscript of a particular vertical value. This array is now scanned for points where the slope between two vertical values is greater than variable TANMAX in absolute value. The value of TANMAX is to be chosen so that it represents the writing rate limit of the digitizer. For any pair found to have a slope greater than TANMAX, the second value of the pair becomes a gap point. It does so if the first value of the pair is valid and the trace could not actually move from there to the second point.

c. Steps a and b have probably produced at least a few gap points. These are filled by linearly interpolating between the vertical values that remain.

d. The gap points at either end cannot be interpolated. Except for the baseline, extrapolation does not seem suitable. The gaps at the ends of the calibration and data traces are set so that they are

---

<sup>2</sup>Robert Puttcamp, SKYNET Applications Software Package (SASP) Programmers' Handbook, Version 2, Harry Diamond Laboratories TR-1844 (February 1978).

zero after steps e and f. The baseline is extended by averaging the five valid vertical values at the end and using that average for the gap points at that end.

e. The baseline trace is subtracted, point for point, from the calibration trace and the data trace.

f. The value 256 is subtracted from each vertical value to place zero at the center line.

The following subroutines are used:

<u>Step</u>	<u>Subroutine</u>	<u>Type of Trace</u>
a	BASATC	Baseline
	CALATC	Calibration
	DATATC	Data
b	ELIMPT	Data
	QKELIM	Baseline, calibration
c	INTRP	Data
	QKLNTP	Baseline, calibration
d,e,f	EXTRAP	Baseline
	ISBASE	Calibration
	SUBASE	Data

#### 4.2 Data Evaluation and Characterization

During the data-taking sequence, four types of data are available for evaluation: baseline, calibration, raw, and normalized. The programs used to plot these data are CHKBAS, CHKCAL, RAWCPY, and MARK, respectively. The STARS operator is responsible for reviewing the data for quality and completeness. He is not typically an experimenter and, therefore, cannot be expected to judge the data in the context of the experiment. Instead, he should comment on data deficiencies or excesses relating to hardware misuses, limitations, or malfunctions.

Data are characterized in task MARK (sect. 4.2.4). There, key parameters of interest to the experimenter are marked with the graphic cursor by the STARS operator.

#### 4.2.1 Baseline Data

Prior to an actual experiment (simulator shot), baseline data must be acquired from each digitizer so that the proper intensity setting for each digitizer may be determined. Through use of the TV monitor, the baseline may be observed while the position and intensity are adjusted. Then program BASE is run to digitize, acquire, normalize, and store the baseline data on a disc. If any intensity settings were incorrect, an appropriate error message appears on the terminal.

Typical error messages are "BLM," bloom too large or too small, and "SHF," baseline trace position not as requested. "BLM" indicates that the intensity setting for that digitizer was either too high or too low. This setting can be corrected by adjusting the intensity and rerunning program BASE. "SHF" is usually caused by an incorrect baseline position. Repositioning the baseline and rerunning program BASE corrects this error. (App A lists and discusses error codes.)

Baseline data are used primarily for two purposes: (1) The baseline waveform is subtracted from the calibration and experiment data waveforms to reform any nonuniformities, such as baseline tilt. (2) Monitoring the baseline is a quick method of determining if the digitizers are functioning properly and the desired knob settings are correct. The turnaround time for program BASE is typically less than 2 min.

If several attempts at acquiring a baseline on one or more digitizers do not produce the desired result, two separate programs are available as diagnostics. Program CHKBAS can be used to graph the baselines on the console display, and program BLMCHK can be used to print the bloom parameter for each digitizer on the operator's console. Once error flags are noted, they must be reset by running program CLEAR before new baselines can be acquired from the suspect digitizers. If a digitizer or channel error flag is not reset, no data from that digitizer or channel are processed or output in any subsequent routine in the current production run. If only one or a few baselines need to be rerun, the other digitizers can be flagged out by running program FLAG and disabling them. Program CLEAR must be run after program BASE to re-enable these digitizers.

#### 4.2.2 Calibration Data

A calibration waveform is obtained from each optical transmitter calibration generator. The calibration pulse amplitude as recorded at the digitizer is a measure of the optical data link gain. In program CHKCAL, calibration waveforms for each channel are plotted on the R4010 terminal. The peak positive or negative value of the calibration pulse is marked by the operator using the graphic cursor on

the R4010 terminal. Then the zero value of the calibration pulse is marked. The data link gain is calculated and stored in disc file MAIN.DAT.

#### 4.2.3 Raw Data

Programs DATAPC, CONDAT, and MORDAT are used to acquire the raw data from the digitizers and store the data on a disc. Raw data are the multivalued functions read directly from the digitizer target (fig. 12). Although they are not useful for mathematical machine computations in this format, the raw data are very useful to the experimenter for a quick look and to the terminal operator in determining if the digitizer writing rate has been exceeded. Program RAWCPY plots the raw data for each digitizer and permits entry of an operator comment prior to producing a hard copy. These raw data are made available to the experimenter within minutes after the shot. At this stage in data evaluation, an experienced operator can predict the success of the normalization in upcoming programs.

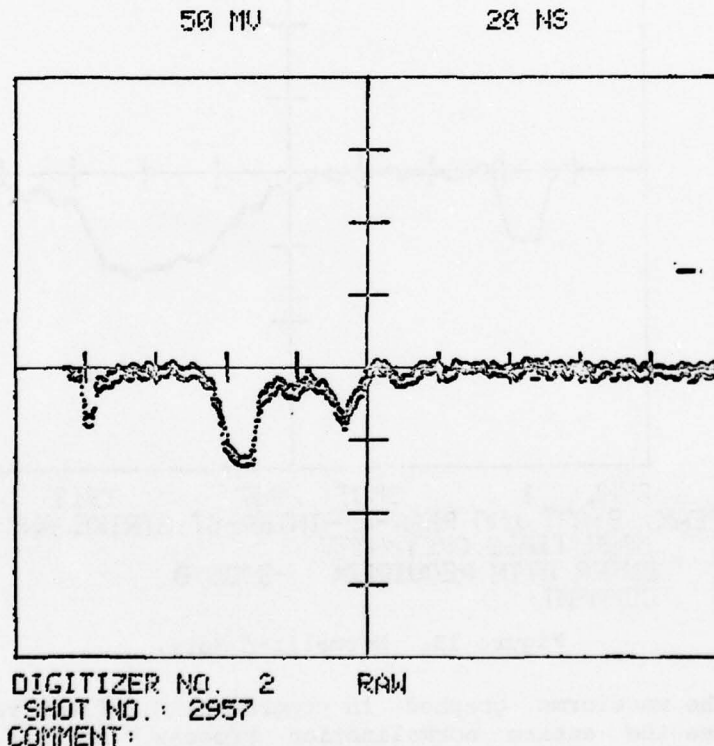


Figure 12. Raw data, sample 1.

#### 4.2.4 Normalized Data

During data taking, the first opportunity to view the normalized data occurs in program MARK. With the normalized, unscaled data graphed on the terminal, the operator evaluates and marks the data waveforms (fig. 13). First, a preprogrammed message is printed instructing the operator to sequence the graphic cursor to three points of interest on the data plot. The time values of these points are then transferred to the main data file on the disc. (Because the data are stored as amplitudes in a 512-element array comprising equal time intervals, knowing the time or array index is sufficient for determining the corresponding amplitude.) Lastly, the prompt appears for a comment to be entered. Comments may be up to 24 characters long.

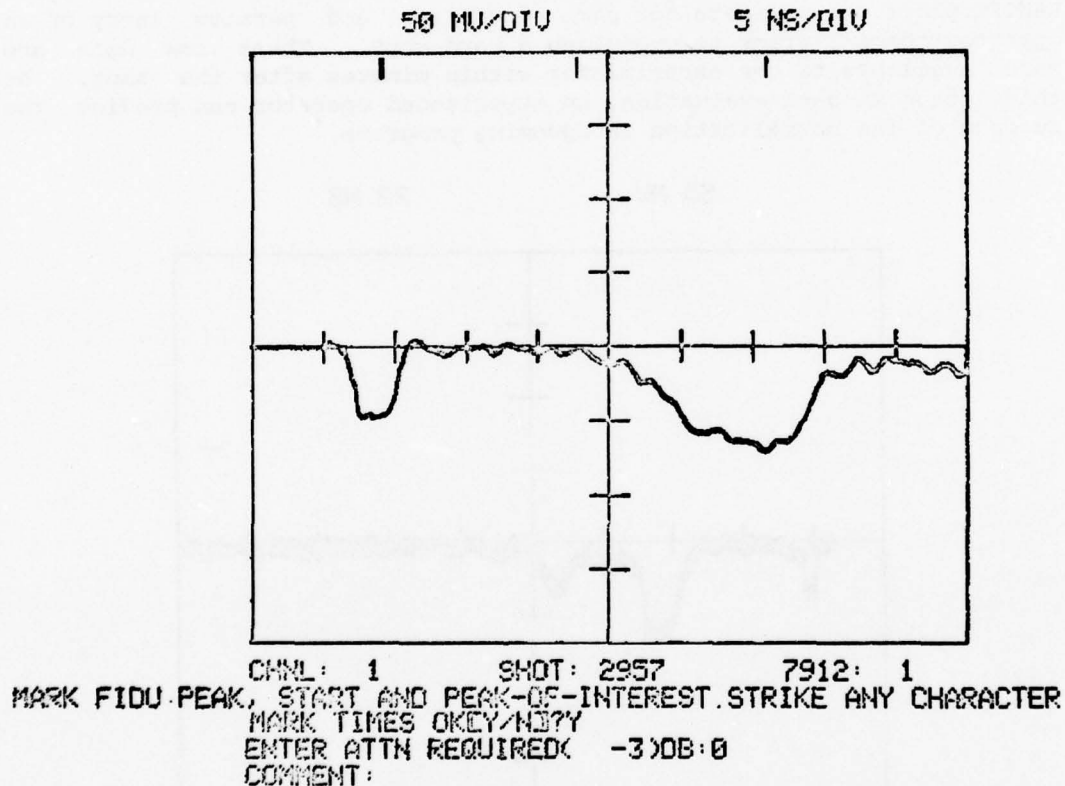


Figure 13. Normalized data.

The waveforms graphed in program MARK reflect raw data that have undergone the entire normalization process. At this stage, hard copies of the raw data have already been printed and can be directly compared one to one with their normalized counterparts in program MARK. Copies of the unscaled, normalized data are not generally retained; however, they could be reproduced later.

### 4.3 Typical Production Run

On each shot, the experimenter is provided with a system configuration printout (fig. 14), which gives the date and time of the event, a synopsis of equipment used on the event, and the output of the two calorimeters associated with the STARS.

SYSTEM CONFIGURATION

TRIGGER DATE AND TIME: 76 8 12 11 46 2

SHOT NO: 2957 2957

NO. OF 7912'S, BOARDS AND CHNLS: 4 2 2

VEN'S & NO. OF BLOCKS: 6 14 22 54 70 71 71

SIM. CODE & DATA: 02 1 0.4220000E 01 0.9999999E-04 0.0000000E 00 0.0000000E 00

COMMENT:

CALORIMETER DATA: 0.1302019E-04 74 0.2021739E-03 100 0

7912	ADDRESS	CHNL	FLAG	SWEEP	VERT	BLOOM	OFFSET
1	164040	1	0	5	50	6	0
2	164050	1	0	20	50	8	0
3	164060	2	0	5	50	4	0
4	164070	2	0	20	50	5	0

CHANL	FLAG	TEST POINT	ATTN	PROBE	BIAS	PRE-X DLY	DEVICE	1ST	2ND	3RD
	XMITTER	X-ATTN	POST-R DLY	POST-R DLY	POST-R GAIN	7912'S				
		GAIN	FACTOR			DELAY				
1	0	ALAU	DML2	I2	0	0.1000000E 01	TRN2			
	16	-12	0.2000000E 01	0.0000000E 00	2	1	2	0		
		0.2527716E 02	0.2316636E-01	0.2203000E 03						
2	0	AGAU	DML1	I1	0	0.1000000E 01	TRN3			
	20	-20	0.2000000E 01	0.0000000E 00	2	3	4	0		
		0.3165972E 02	0.6001152E-01	0.2233000E 03						

Figure 14. System configuration printout.

Within 5 min after the event, the experimenter is given a copy of the raw data for each digitizer (fig. 12, 15). He can quickly estimate the magnitude by dividing the number of digitizer volts by the "factor" given in figure 14. The raw data are annotated with the shot number and the digitizer parameters.


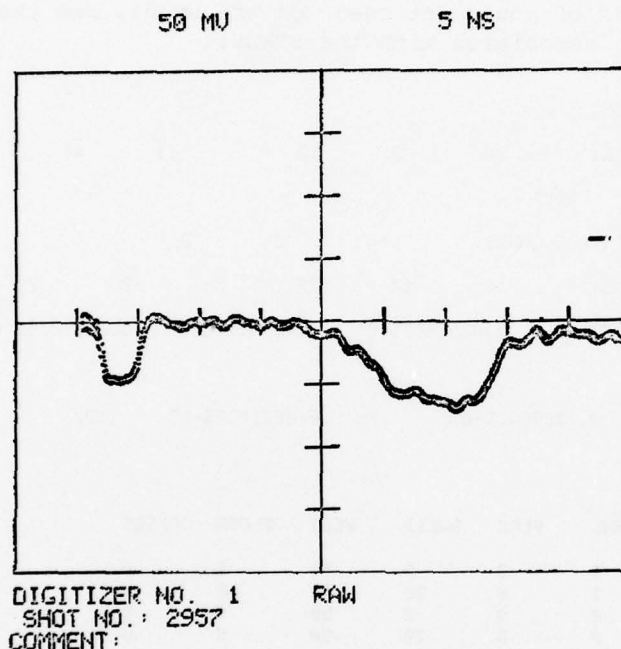


Figure 15. Raw data, sample 2.

Within 45 min, the experimenter is given a copy of the scaled data (fig. 16) in the units desired by him. Accurate timing between channels is possible by referencing the marker pulse (fiducial) at the left edge of all data traces. The STARS operator marks each trace with the location of the fiducial marker, the start time of the data, and the time of the peak of interest. Vertical timing tick marks indicate the exact locations of these parameters on the scaled data plot. Annotation of several key parameters, such as the coded test point name, also is provided for each plot.

Special output routines permit printing of the raw data, the calorimeter sample values, the baseline data, etc. Routines have been created also to allow for printing the contents of each of the calibration files; for example, figure 17 was generated by program PRB. Program PRB lists all the probes currently in the PRBCAL.DAT file and pertinent parameters for each probe.

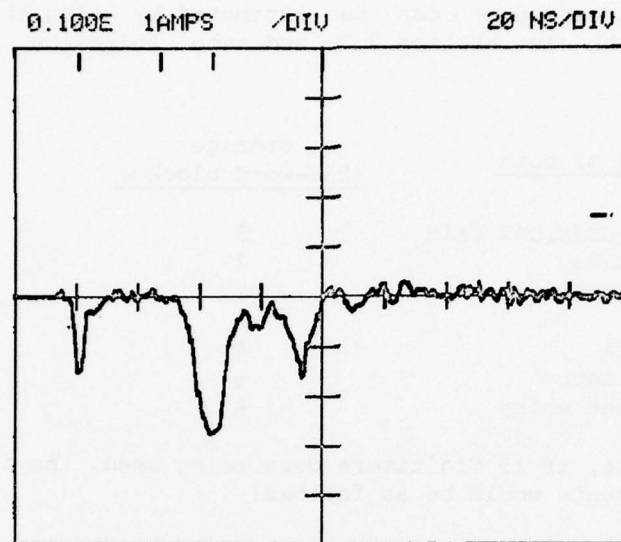

CHNL: 1 SHOT: 2957  
 TEST PT: ALAU 12/ 8/76 11:46: 2  
 PEAK= -27.316E-01AMPS  
 TIME TO PEAK= 17.187E 00NS  
 CALORIMETER TP3: 1.30202E-05CAL/CMMX2  
 CALORIMETER 2B: 2.02174E-04CAL/CMMX2

Figure 16. Scaled data.

ID	XFER(DB)	DLY(NS)	FLAG	BLN	UNITS	MT	PARAMETR	MODEL	MFR S/N	CONNECTOR
B501	-140.03	2.63	1	1	GAUSS	BD	0.10E-02	CMLX5AR	1	2 SMA 50
B506	-140.03	20.00	1	1	GAUSS	BD	0.10E-02	CMLX5AR	6	2 SMA 50
B310	-146.06	10.75	1	1	GAUSS	BD	0.50E-03	CMLX3AR	10	2 SMA 50
B301	-146.06	2.45	1	1	GAUSS	BD	0.50E-03	CMLX3AR	1	2 SMA 50
B507	-140.03	20.00	1	1	GAUSS	BD	0.10E-02	CMLX5AR	7	2 SMA 50
EAU	-5.93	0.00	0	2	V/M	E	0.50E-02	PMD-1	13	1 SMA 50
I1	0.00	0.30	0	2	AMPS	I	0.50E 00	CT-2	1	1 SMA 50
I2	0.00	0.30	0	2	AMPS	I	0.50E 00	CT-2	2	1 SMA 50
B308	-146.06	10.88	1	1	GAUSS	BD	0.50E-03	CMLX3AR	8	2 SMA 50
D309	-146.06	9.25	1	1	GAUSS	BD	0.50E-03	CMLX3AR	7	2 SMA 50
EAG	-5.93	19.00	0	2	V/M	E	0.50E-02	PMD-1	17	1 SMA 50
EGLY	-5.93	19.50	0	2	V/M	E	0.50E-02	PMD-1	16	1 SMA 50
OPEN	0.00	0.00	0	1	VOLTS	V	0.00E 00			

TOTAL NO. OF PROBES: 13 00:35:57 11-AUG-76

Figure 17. Probe file listing.

#### 4.4 Storage and Retrieval

Storage requirements can be estimated by using the format of the shot file given in section 2.2 and the following file entry information:

<u>Type of data</u>	<u>Storage (256-word blocks)</u>
Copy of SYSCON.DAT file	5
Baseline data	2
Calibration data	2
Raw data	8
Scaled data	4
Parameter entry	$\frac{1}{4}$
User comment entry	$\frac{1}{4}$

For example, if 15 digitizers were being used, the file size or disc space requirements would be as follows:

$$\begin{aligned}5 + 2n + 2n + 8n + 4n + n/4 + n/4 &= 5 + 16.5n \\&= 5 + 16.5 \times 15 \\&= 252.5 \text{ blocks .}\end{aligned}$$

This must be rounded up to a whole block value; hence the file would actually use 253 disc blocks.

Further, assume that 10 simulator shots (10 SHOT.DAT files) are to be recorded:  $10 \times 253 = 2530$  blocks. The disc pack used to store these data contains 1,228,800 words. Therefore, a disc pack contains 4800 blocks. For a typical test period of several weeks, some other form of off-line storage is indicated. The STARS uses nine-track tape compatible with that used in industry.

At the end of each day's testing, the data disc contents are transferred to a magnetic tape. This tape then becomes the master data tape. This master data tape is then copied to provide a backup. A test program typically lasts from 2 to 3 weeks. The data from this program should fit on a single reel of tape. Backup copies of all master data file tapes are retained. These tapes then comprise the STARS data bank.

Certain features have been incorporated into the shot file structure to facilitate searching. In the SYSCON header entry of the shot file (fig. 3), virtual block number (VBN) entries point to the

first entry of each type shown. For example, on shot No. 2937 for digitizer No. 1, the sequence for obtaining the raw data would be as follows:

a. The operator locates the master file tape containing shot No. 2937 and mounts this tape on the nine-track tape unit.

b. The operator uses the appropriate program to search for the proper file. For example, the shot file may be "S293702.DAT." The "02" is the simulator code for the OWL II simulator. The "S" and the "DAT" are always the same.

c. The operator transfers the data to the disc.

d. The SYSCON header in the S293702.DAT file contains an entry entitled "VBN, RAW" which is an integer representing the block number within the S293702.DAT file of the raw data of digitizer No. 1. Reading forward eight blocks transfers the data of interest. (Raw data are always eight blocks in length.)

Another example is the test point name contained in the SYSCON channel entry of the shot file. The test point name permits nine ASCII characters that can be coded into representing some important parameters of the test. For example, three digits can be used to code 1000 different location or test point descriptions, two digits can denote the measurement type, and two digits can denote the simulator orientation and type. This type of coding facilitates search algorithms.

## 5. BRIEF MODULE DESCRIPTIONS

This section presents a short synopsis of each program and subroutine in SASP. Detailed descriptions of each module can be found in the SASP Programmers' Handbook.<sup>2</sup> Indirect command file TESTSEQ is discussed in section 3.4.1. Wherever possible, the subroutines and macros have been grouped logically into object and macro libraries resident on the SASP disc (app E).

### Indirect command file

### Function

TESTSEQ	Control data acquisition and processing program execution.
---------	--

### Task

### Function

ATN	Generate, modify, and print contents of ATNCAL.DAT file.
-----	--

<sup>2</sup>Robert Puttcamp, SKYNET Applications Software Package (SASP) Programmers' Handbook, Version 2, Harry Diamond Laboratories TR-1844 (February 1978).

<u>Task</u>	<u>Function</u>
BASE	Acquire, normalize, and store baseline for each enabled digitizer.
BLKCNT	Print number of blocks required for shot file (that is, MAIN.DAT and Snnnnxx.DAT).
BLMCHK	List digitizer entries in SYSCON COMMON area in core on operator's console.
CAL	Acquire, normalize, and store calibration trace for each enabled digitizer.
CHKBAS	Plot baseline data.
CHKCAL	Plot and evaluate calibration waveforms, compute optical link gain, and input comments.
CLEAR	Reset digitizer and channel flags to values stored in disc file FLAG.DAT.
CONDAT	(See DATAPC.)
CONGEN	Generate SYSCON.DAT file and initialize digitizer addresses.
DATAPC, CONDAT, MORDAT	Acquire and store raw data trace for each enabled digitizer; acquire, analyze, and print calorimeter data; these programs are always run as chain.
ERASE	Erase screen of terminal or give eight line feeds on LA36 DEC Teletype terminal.
FIXNPT	Allow partial recovery from overflow condition on data trace; fix data point count of raw data at 2000; on overflow, store first 2000 points read.
FLAG	List digitizer and channel error flag and allow enabling or disabling.
FPLOT	Plot final scaled data produced by program VCOPY
GETCON	Read file SYSCON.DAT into COMMON area SYSCON.
KNOB	Switch all enabled digitizers to knob setting mode (TV and graticule).
LOG	Input log information to SYSCON.DAT file.

<u>Task</u>	<u>Function</u>
MARK	Plot and characterize normalized data waveforms and input comments.
MORDAT	(See DATAPC.)
OFFLIN	Test each enabled digitizer for off-line condition.
OUTFIL	Create MAIN.DAT file.
PLOTCL	Dummy FORTRAN program to call PLTCLR to get around some system problems with DEC FORTRAN.
PLTCLR	Plot calorimeter data on R4010 terminal and make hard copy of them.
PRB	Generate, modify, and print contents of PRBCAL.DAT.
PREDAT	Prepare enabled digitizers for data taking (required precursor for DATAPC, CONDAT, and MORDAT).
PRTNAM	Print name of shot file to which SYSCON COMMON area points.
RAWCPY	Plot, annotate, and comment on raw data waveforms.
RECCAL	Reproduce calibration traces and operator's marks on R4010 terminal and make hard copy of plots.
SAVCON	Write file SYSCON.DAT from COMMON area SYSCON.
SCACPY	Plot, annotate, and rescale data waveforms using R4010 terminal and hard-copy device.
SCALE	Normalize and scale raw data traces, calculate overall scale factor and data link gain, and store results.
SCAVP	Plot, annotate, and rescale data waveforms using Versatec plot package and M1110 printer plotter.
SETFLG	Generate master flag setting file FLAG.DAT.
SHOTNO	Input simulator code and shot number to SYSCON COMMON area.
SNAPIT	Dump SYSCON common partition to line printer in standard RSX-11M format.
SUMERR	Print digitizer and channel error flags.

<u>Task</u>	<u>Function</u>
SYSCHK	Check validity of certain entries in SYSCON COMMON area.
SYSDMP	Quickly dump contents of SYSCON COMMON area.
TABIN	Input simulator information to SYSCON.DAT file, and rename MAIN.DAT file to current shot file name.
TABLES	Generate three tables of user-oriented tabular output data.
USRCOM	Input user comments to shot file.
VCOPY	Generate disc file containing final scaled data graphs.
XCH	Generate, modify, and print contents of XHCAL.DAT file.

<u>Subroutine</u>	<u>Function</u>
ABSPKV	Compute absolute peak value of ATC array.
ANOCLE	Annotate plots of calorimeter data.
ANOMRK	Annotate MARK (unscaled, normalized data).
ANOSCA	Annotate parameters on scaled data.
ANOTS	Annotate Versatec plot package final data graphs.
ANSWER	Input and decode "yes" or "no" answer from console.
ASCBIN	Convert ASCII input to binary integers.
ASDCD	Transfer selected ASCII field from input buffer.
ASTSC	Set and clear asynchronous system trap (AST) for RAWPLT program; recognize ESC key on TTL:.
ATNREF	Handle reference to ATNCAL.DAT file.
ATNSCH	Read ATNCAL.DAT file for entry of selected device.
ATTEN	Set and verify remote attenuator settings, produce baseline, and calibrate triggers.
BASATC	Perform ATC for baseline trace, including bloom calculation.
BUSY	Wait for clearing of busy status of digitizer.

<u>Subroutine</u>	<u>Function</u>
CALATC	Perform ATC for calibration trace.
CALCSF	Calculate overall scale factor and data link gain.
CALIB	Control transmitters during acquisition of calibration trace.
CALOR	Process data from calorimeters.
CALPLT	Actually plot calorimeter data on R4010 terminal.
CALSPC (CALIB)	Control transmitters during acquisition of calibration trace, but do not flag out either digitizers or channels for errors.
CHACNT	Count number of nonblank ASCII characters in integer array.
CHKTRG	Check each enabled digitizer for trigger interrupt.
CHNERR	Flag channel errors.
CLRCBD	Clear digitizer controller boards' status words, thereby disabling interrupts.
CMDDIG	Send selected command to all enabled digitizers and confirm its execution.
DATATC	Perform ATC for data trace.
DATFIO	Handle input and output to MAIN.DAT and shot files.
DCDBIT	Generate table of codes corresponding to bits set in given word.
DIGEFR	Flag digitizer errors.
DISABL	Set disabled flag in digitizer or channel error word.
DSKERR	(AST routine) Handle disc errors.
ELIMPT	Eliminate noise points from data trace of ATC.
ERRNUM	Return bit number of left-most bit set in given word.
EXTRAP	Extrapolate ends of baseline trace, and subtract 256 from entire baseline trace.
FIELD	Set pointer to selected field in input buffer, and count characters in field.

<u>Subroutine</u>	<u>Function</u>
FLDCD	Decode floating point number from selected field in input buffer.
FLGPTS	Flag raw data points known to be incorrect.
FLOTE	Convert integer to floating point number.
FLPOP	(Multientry point subroutine) Handle floating point arithmetic.
FSCIN	Read copy of SYSCON from shot file into SYSCON COMMON.
FXDCD	Decode integer from selected field in input buffer.
GETBAS	Read, reduce, and store baseline trace for all enabled digitizers.
GETCAL	Read, reduce, and store calibration trace for all enabled digitizers.
GETDAT	Read and store raw data trace for all enabled digitizers.
GRAATC	Graph normalized floating point data.
GRABAS	Graph baseline data.
GRACAL	Graph calibration data.
GRAPHV	Plot scaled data for final data plots from Versatec plot package.
GRARAW	Plot raw data.
GRASCA	Plot scaled data.
GRAT	Plot graticules on terminal in three sizes.
IASCII	Convert integer of less than five digits to ASCII.
IFICKS	Convert floating point number to integer (truncate).
INITLZ	Read and store data from calorimeters.
INPCSR	Input data characterization parameters with graphic cursor.
INPUT	Input line from logical unit number (LUN) 5 to input buffer.

<u>Subroutine</u>	<u>Function</u>
INT	Integrate scaled data by trapezoidal integration method.
INTABL	Enable digitizer interrupts.
INTDIG	Handle digitizer interrupts.
INTRP	Fill gaps in data trace by linear interpolation.
ISBASE	Subtract baseline trace and 256 from calibration trace.
LSTERR	List channel or digitizer errors.
MARKIT	Put cursor marks generated by program MARK on final data graphs for Versatec plot package.
MINMAX	Search integer array for minimum and maximum values.
MRKCUR	Print cursor marks generated by subroutine INPCSR on scaled and unscaled data plots.
NSCALE	Adjust floating point number to value between 1 and 10, and calculate proper exponent of 10 to keep same value.
PAGE	Erase screen at terminal or give eight line feeds on LA36 Teletype terminal.
PLTSF	Annotate digitizer scale factors on baseline, calibration, ATC (normalized), or raw data plots.
PRBREF	Handle references to PRBCAL.DAT file.
PRBSCH	Read PRBCAL.DAT file for entry of selected probe.
PRICLR	Print raw calorimeter data.
QKELIM	Eliminate noise points from ATC of baseline or ATC of calibration trace.
QKLNTP	Fill gaps in baseline or calibration trace by linear interpolation.
RAWPLT	Condition raw data array for plotting by subroutine GRARAW.
RDIG	Read raw data from digitizer.
RDUKRO	Reduce knob readout information in raw data and compare it with desired settings.

<u>Subroutine</u>	<u>Function</u>
SAVCLR	Write calorimeter data into file on DK1: having same name as shot file, but starting with C instead of S.
SCADAT	Normalize and scale raw data traces, and store results.
SCGEN	Create SYSCON.DAT file, and write initial information.
SCIO	Read and rewrite SYSCON.DAT file into and from SYSCON COMMON area in core.
SCMAIN	Copy SYSCON.DAT file into MAIN.DAT file, and rename MAIN.DAT file to current shot file name.
SETVEC	Initialize digitizer interrupt vectors.
STCDAT	Transfer trigger date and time from TIME variable to SYSCON COMMON area in core.
SUBASE	Subtract baseline trace and 256 from data trace.
UCOUT	Write user comments into shot file.
VGRAT	Produce graph frame, cross bars, and tick marks for graphs produced with Versatec package.
XCHDRI	Read header records of XCHCAL.DAT file.
XCHREF	Handle references to XCHCAL.DAT file.
XCHSCH	Read XCHCAL.DAT for entry of selected channel.
XMIT	Control transmitters during acquisition of data trace.
XMTSPC	Control transmitters during acquisition of data trace; report, but do not flag out, attenuator verification errors.
.RSTR	Restore registers from stack.
.SAVE	Save registers on stack.
<u>Macro</u>	<u>Function</u>
RDSECT	Establish PSECT (program section) for RAWDAT COMMON area
SCSECT	Establish PSECT for SYSCON COMMON area.
SYMDEF	Define commonly used symbols.

#### ABBREVIATIONS

ASCII	American Standard Code for Information Interchange
ATC	Average to center of trace
CR	Carriage return
CRO:	Card reader unit No. 0
CT0:	Cassette tape unit No. 0
CT1:	Cassette tape unit No. 1
DEC	Digital Equipment Corp.
DK0:	Disc unit No. 0, system disc
DK1:	Disc unit No. 1
LA36	DEC Teletype terminal (TT1:)
LP0:	Versatec M1110 line printer, unit No. 0
LUN	Logical unit number
MCR	Monitor console routine
MT0:	Magnetic tape unit No. 0
M1110	Printer plotter
OWL II	Oil Water Line II, SGEMP simulator
RK11	DEC type number for disc drives on PDP 11/40
RSX-11M	DEC operating system
R632	Tektronix television monitor
R4010	Tektronix main console graphic display terminal (TT0:)
R7912	Tektronix transient digitizer, primary data recording device
SASP	SKYNET Applications Software Package

ABBREVIATIONS (Cont'd)

SGEMP	Systems generated electromagnetic pulse
STARS	SGEMP Transportable Automated Recording System
SY:	System disc, usually DK0:
TI:	Operator's current console, either TT0: or TT1:
TV	Television
UFD	User file directory
UIC	User identification code
VBN	Virtual block number
XRD	X-ray detector

# APPENDIX A.--DIGITIZER AND CHANNEL ERROR CODES

## A-1. SUMMARY

Both the digitizer and channel entries in the SYSCON.DAT file contain integer words for error flags, each bit of which is coded as shown below.

### Digitizer Error Words

<u>Bit No.</u>	<u>Octal equivalent</u>	<u>Code</u>	<u>Meaning</u>
0	000001	OFF	Off line
1	000002	DIS	Disabled
2	000004	OVF	Data overflow
3	000010	UNF	Data underflow
4	000020	BZY	Busy status
5	000040	CMD	Command execution error
6	000100	TRG	Failure to trigger
7	000200	VFY	Reading verification error
8	000400	KRO	Bad (incorrect) knob readout
9	001000	BLM	Bloom too large or too small
10	002000	SHF	Baseline trace position not as requested
11	004000	(Unassigned)	
12	010000	(Unassigned)	
13	020000	NOB	No baseline trace
14	040000	NOC	No calibration trace
15	100000	NOD	No raw data trace

## APPENDIX A

Channel Error Words

<u>Bit No.</u>	<u>Octal equivalent</u>	<u>Code</u>	<u>Meaning</u>
0	000001	OFF	Offline
1	000002	DIS	Disablement
2	000004	ATN	Attenuator setting verification error
3	000010	CAL	Bad calibration trace
4 to 5	(Unassigned)		

## A-2. CAUSES AND RECOVERY

Following is a list of the digitizer and channel error codes, the conditions which cause the error to be generated, and the procedure to correct the faulty condition.

<u>Digitizer error code</u>	<u>Cause</u>	<u>Recovery</u>
0000	1. Digitizer is off or physically disconnected.	1. Turn on or reconnect digitizer.
0001	1. Digitizer is intentionally disabled (with program FLAG). 2. Error is on associated channel.	1. Reenable digitizer with program FLAG. 2. Reenable digitizer with program FLAG.
0002	1. Too many points are stored in digitizer memory.	1. Decrease intensity or increase sweep rate (faster sweep speed).
0003	1. Too few points are stored in digitizer memory.	1. Ascertain proper trigger available, increase intensity, or both.
0004	1. Busy status exists.	1. Allow digitizer to warm up. Manually switch to digital and back to non-store on digitizer front panel.

## APPENDIX A

<u>Digitizer error code</u>	<u>Cause</u>	<u>Recovery</u>
0005	<ol style="list-style-type: none"> <li>1. Command is not executed.</li> <li>2. Digitizer is not in proper status (such as single shot mode).</li> </ol>	<ol style="list-style-type: none"> <li>1. Check digitizer front panel and plug-in switches for nonstore, single sweep, and external trigger.</li> </ol>
0006	<ol style="list-style-type: none"> <li>1. Trigger cable is disconnected.</li> <li>2. Trigger level and slope are improper.</li> </ol>	<ol style="list-style-type: none"> <li>1. Reconnect cable.</li> <li>2. Readjust level and slope.</li> </ol>
0007	<ol style="list-style-type: none"> <li>1. On variable TRYMAX separate attempts, second reading of digitizer memory had more than variable MISMAX words different from first reading.</li> </ol>	<ol style="list-style-type: none"> <li>1. Repair digitizer after repeated errors</li> </ol>
0008	<ol style="list-style-type: none"> <li>1. Setting of plug-in does not match that requested.</li> <li>2. American Standard Code for Information Interchange (ASCII) readout is incorrect.</li> </ol>	<ol style="list-style-type: none"> <li>1. Change request in SYSCON.DAT file with program LOG or change plug-in setting.</li> <li>2. Repair ASCII board. Temporarily reenable it and continue as though no error occurred. Ascertain that plug-ins are set as desired.</li> </ol>
0009	<ol style="list-style-type: none"> <li>1. Intensity is too high or too low.</li> </ol>	<ol style="list-style-type: none"> <li>1. Adjust intensity and rerun program BASE.</li> </ol>
0010	<ol style="list-style-type: none"> <li>1. Baseline position is incorrect.</li> <li>2. Noise rides on baseline.</li> <li>3. Noise point in digitizer is not eliminated by normalization.</li> </ol>	<ol style="list-style-type: none"> <li>1. Reposition baseline.</li> <li>2. Check input signal.</li> <li>3. Rerun program BASE.</li> </ol>

# APPENDIX A

## Digitizer error code

## Cause

## Recovery

- |      |  |   |
|------|--|---|
| 0013 | 1. There was a previous error.                       | 1. If previous error occurred before acquisition and reduction of baseline trace, there is truly no baseline; therefore, rerun program BASE after correcting error.<br><br>2. If previous error was 0008, 0009, or 0010, baseline data are stored and can be checked. Correct error and rerun program BASE. |
| 0014 | 1. There was a previous error during program CAL.    | 1. Follow directions for error 0013, using task CAL.  |
| 0015 | 1. There was a previous error during program MORDAT. | 1. Follow directions for error 0013, using program MORDAT.  |

## Channel error code

## Cause

## Recovery

- |      |  |   |
|------|--|---|
| 0000 | 1. Transmitter is not connected.   | 1. Connect transmitter.   |
| 0001 | 1. Channel was intentionally disabled.<br><br>2. Associated digitizer had error. | 1. Reenable channel.  |
| 0002 | 1. Transmitter attenuator did not respond properly.                              | 1. Manually override attenuator..<br><br>2. Repair attenuator.                            |
| 0003 | 1. Operator considered calibration trace to be wrong.                            | 1. Obtain good calibration trace with program CAL, and recheck trace with program CHKCAL. |

A-3. GENERAL

An error condition automatically prevents the unit and its associated data from being accessed further. All programs examine the error word. Before any form of recovery is attempted or data are accessed, the operator should use program CLEAR to reenable the unit.

If he desires to eliminate a digitizer from consideration, the operator should use program FLAG to disable the digitizer, thereby creating error condition 0001. To permanently disable a digitizer, use program SETFLG.

An error condition generates a disable error (0001) also on the associated units. The exception occurs when only one of the digitizers connected to a channel is in error. Under this condition, the channel is not disabled.

## APPENDIX B.--DIGITIZER COMMANDS AND STATUS

A fairly comprehensive description of programming for the Tektronix R7912 transient digitizer can be found in the appropriate Tektronix manual;<sup>1</sup> only the basics are described here.

These Digital Equipment Corp. (DEC) UNIBUS<sup>2</sup> addresses are used in programming:

16X000 <sub>8</sub>	Digitizer controller status word for card No. (X-3)
16X0D2 <sub>8</sub>	Digitizer command register for digitizer defined as having address 16X0D0 <sub>8</sub>
16X0D4 <sub>8</sub>	Digitizer data and status register for digitizer defined as having address 16X0D0 <sub>8</sub>

The controller status word content is defined as follows:

<u>Bit</u>	<u>Definition</u>
0	Interrupt pending
3,4,5	Device number ("D" of address) of interrupting device
6	Interrupt enable

The digitizer status word content is defined as follows:

<u>Bit</u>	<u>Definition</u>
0	Identity as R7912 digitizer
1	Remote channel mode
2	Memory locked (write protect)
3	Graticule enabled
4	Single sweep armed
5	Single sweep mode

<sup>1</sup>WDI TEK BASIC Software User's Manual, Tektronix, Inc., Beaverton, OR (1974).

<sup>2</sup>DEC PDP 11 Processor Handbook, 1975-1976, Digital Equipment Corp., Maynard, MA (1975).

## APPENDIX B

<u>Bit</u>	<u>Definition</u>
6	Channel 1 selected (if in remote)
7	Digital mode
8 to 12	(Not used)
13	Busy

The following are the digitizer commands where "X=0" equates to "clear," and "X=1" equates to "enable":

<u>Command</u>	<u>Function</u>
12000 <sub>8</sub>	Load memory, repetitive sweep only.
20X2 <sub>8</sub>	Inhibit knob readout output.
20X3 <sub>8</sub>	Select channel No. 1.
20X4 <sub>8</sub>	Select dot graticule.
20X5 <sub>8</sub>	Select digital mode.
20X6 <sub>8</sub>	Select remote channel mode.
20X7 <sub>8</sub>	Lock memory.
2020 <sub>8</sub>	Read status (status placed in status register).
2040 <sub>8</sub>	Arm single sweep.
2100 <sub>8</sub>	Read data (data word placed in data register as previous word moved out).
4000 <sub>8</sub>	Initialize registers.

For direct diagnostic control, to enter commands at the front panel switches of the central processing unit, the operator deposits the front panel switch contents to the command register address. To read the status, the operator examines the status register. The DEC PDP 11/40 Processor Handbook<sup>2</sup> describes the console operation. With the memory management option in place, the direct addresses are 76X0D0<sub>8</sub> instead of 16X0D0<sub>8</sub> although the virtual address for privileged tasks remains as 16X0D0<sub>8</sub>.

<sup>2</sup>DEC PDP 11 Processor Handbook, 1975-1976, Digital Equipment Corp., Maynard, MA (1975).

## APPENDIX B

More complete diagnostic capabilities exist in the Waveform Digitizing Instrument (WDI) TEK BASIC operating system loadable from a cassette or disc. Both the cassette and the disc are available at the Systems Generated Electromagnetic Pulse Transportable Automated Recording System (STARS) facility to aid in troubleshooting R7912 hardware problems.

#### APPENDIX C.--TESTSEQ OPERATIONS CONTROL FILE

The operator of the Systems Generated Electromagnetic Pulse Transportable Automated Recording System (STARS) can control the data-taking sequence and the execution of all programs by an indirect command file, TESTSEQ.CMD. This file is composed of all the programs used to acquire, process, and store the data. When the operator types "@DK1:TESTSEQ" and then a carriage return at the console terminal, the complete directions for taking data are printed at the console display. The instructions lead the operator through the entire sequence and permit interactive control of program execution.

Following is the TESTSEQ file. Reader familiarity with the command language of the monitor console routine of the Digital Equipment Corp. RSX-11M operating system is a prerequisite to understanding the sample.

```
RUN ERASE
;THE FOLLOWING IS AN OUTLINE OF THE DATA TAKING SEQUENCE
;
;EVALUATE DISK SPACE AVAILABILITY (BLKCNT, VFY)
;CREATE MAIN.DAT (OUTFIL)
;RESET 7912/CHANNEL ERROR FLAGS (FLAG)
;CHECK FOR 7912'S OFF-LINE (OFFLIN)
;ENTER LOG INFO (LOG)
;CHECK LOG ENTRIES (SYSCHK)
;CORRECT ANY BAD ENTRIES (LOG, PRB, ATN)
;SET KNOBS (OPTIONAL TV MODE (KNOB))
;TAKE BASELINE (BASE)
;TAKE CAL (CAL)
;EVALUATE CAL (CHKCAL)
;ENTER SHOT NO. AND SIMULATOR CODE (SHOTNO)
;PREPARE 7912'S FOR DATA TAKE (PREDAT)
;TAKE DATA (DATA)
;NORMALIZE AND SCALE DATA (SCALE)
;HARD COPY RAW DATA (RAWCPY)
;EVALUATE DATA QUALITY AND SELECT KEY POINTS (MARK)
;INPUT SIMULATOR INFORMATION (TABIN)
;HARD COPY SCALED DATA (SCAPY)
;PRINT PARAMETER TABLE (TABLES)
;PRINT SYSTEM DESCRIPTION TABLES (TABLES)
;INPUT USER COMMENTS (USRCOM)
;PUT SHOT FILE ON DAILY DATA TAPE (PRTNAM, FLX)
;
;TASK NAMES ARE IN PARENTHESES.
```

# APPENDIX C

```

;MAKE A COPY OF THIS PAGE IF ONE IS NOT AT HAND
;PAUSE
RUN ERASE
;TO ELIMINATE SUPPLEMENTARY, EXPLANATORY TEXT,
;ANSWER THE FOLLOWING WITH N.
;ASK DIALOG EXTENDED DIALOGUE
;
.100:      ;DISK SPACE REQUIREMENTS
RUN BLKCN:
;IFNINS...VFY      INS...VFY
.XQT  VFT  TTO: = DK1:/FR
.WAIT
;
.ASK SPACE IS THERE ENOUGH SPACE
.IFT SPACE      .GOTO 200
.IFF DIALOG      .GOTO 110
;DURING THE FOLLOWING PAUSE,
;MAKE USE OF PIP, FLX AND OTHER UTILITIES TO
.110:      ;CLEAN UP THE DATA DISK
.PAUSE
.GOTO 100

.200:      RUN OUTFIL
.IFT DIALOG:MAIN.DAT HAS NOW BEEN CREATED
;
;PUT ATTENUATORS IN AUTOMATIC MODE
.SETF REPCAL
.SETF REPRED

.300:      RUN FLAG
RUN ERASE
.IFT DIALOG: NOW CHECKING FOR 7912'S OFF-LINE
RUN OFFLIN
.IFF DIALOG      .GOTO 310
;IF PRESENT ERROR CONDITIONS ARE NOT SATISFACTORY
;ANSWER THE FOLLOWING WITH N.
.310:      .ASK GOOD  GOOD STATUS
.IFT GOOD      .GOTO 400
.IFF DIALOG      .GOTO 320
;DURING THE FOLLOWING PAUSE,
;RESOLVE THE STATUS OF THE HARDWARE
;AND PREPARE DECISIONS FOR ENABLE/DISABLE.
.320:      .PAUSE
.GOTO 300

.400:      .IFT REPCAL      .GOTO 800
.IFT REPRED      .GOTO 900
RUN LOG

```

## APPENDIX C

```
.500:      RUN SYSCHK
        .IFF DIALOG .GOTO 510
        ;IF ANY ENTRIES WERE BAD,
        ;ANSWER THE FOLLOWING WITH N.
.510:      ;ASK GOOD ALL ENTRIES GOOD
        .IFT GOOD .GOTO 600
        RUN ERASE
        .IFF DIALOG .GOTO 520
        ;DURING THE FOLLOWING PAUSE,
        ;MAKE USE OF LOG, PRB AND ATN
        ;TO REMOVE THE BAD ENTRIES.
        ;UPON RESUME, ENTRIES WILL BE RECHECKED.
.520:      .PAUSE
        .GOTO 500

.600:      RUN ERASE
        .IFF DIALOG .GOTO 610
        ;IF USE OF THE TV MONITORS
        ;FOR KNOB SETTING IS DESIRED,
        ;ANSWER THE FOLLOWING WITH Y.
.610:      .ASK TV TV MODE
        .IFF TV .GOTO 700
        RUN KNOB
        .IFT DIALOG; 7912'S NOW IN TV MODE
        ;SWITCH TO LOWER INTENSITY
        ;THEN USE INTERNAL TRIGGER

.700:      ;SET KNOBS
        .IFT TV; RESET TRIGGER AND INTENSITY
        .IFF DIALOG .GOTO 710
        ;USE THE FOLLOWING PAUSE.
        ;RESUME WILL INITIATE BASELINE SHOT.
.710:      .PAUSE
        RUN BASE
        .SETF REPCAL
        .SETF REPRED
        .IFF DIALOG .GOTO 720
        ;IF UNDESIRABLE ERRORS HAVE OCCURRED,
        ;ANSWER THE FOLLOWING WITH N.
        ;A Y WILL BRING A CAL SHOT IMMEDIATELY
.720:      .ASK GOOD SATISFACTORY RESULTS
        .IFF GOOD .GOTO 300

.800:      RUN CAL
        RUN CHKCAL
        RUN ERASE
        .IFF DIALOG .GOTO 810
        ;IF THE CAL DATA ARE NOT ACCEPTABLE
```

# APPENDIX C

```

;OR UNWANTED ERRORS HAVE OCCURRED,
;ANSWER THE FOLLOWING WITH N.
.800:      .ASK GOOD SATISFACTORY RESULTS
          .IFT GOOD   .GOTO 900
          .IFF DIALOG .GOTO 820
          ;IF CHANGES TO SYSCON (OTHER THAN FLAGS)
          ;NEED NOT BE MADE AND IF A BASELINE
          ;SHOT NEED NOT BE REPEATED,
          ;ANSWER THE FOLLOWING WITH Y.
.820:      .ASK REPCAL REDO CAL ONLY
          .SETF REPRED
          .GOTO 300

.900:      RUN ERASE
          RUN SHOTNO
          RUN ERASE
          RUN PREDAT
          .IFF DIALOG .GOTO 910
          ;7912's NOW READY FOR DATA RUN
          ;IF ANY UNWANTED ERRORS HAVE OCCURRED,
          ;ANSWER THE FOLLOWING WITH N.
.910      .ASK GOOD SATISFACTORY RESULTS
          .IFT GOOD   .GOTO 1000
          .IFF DIALOG .GOTO 920
          ;IF CHANGES TO SYSCON (OTHER THAN FLAGS)
          ;NEED NOT BE MADE AND IF BASELINE AND
          ;CAL NEED NOT BE REPEATED,
          ;ANSWER THE FOLLOWING WITH Y.
.920:      .ASK REPRED REDO PREDAT ONLY
          .SETF REPCAL
          .GOTO 300

.1000:     RUN ERASE
          .IFF DIALOG .GOTO 1010
          ;A SHORT TIME AFTER ANSWERING THE
          ;FOLLOWING QUESTION WITH Y, THE SYSTEM
          ;WILL BE AWAITING A TRIGGER PULSE.
          ;ANSWER WITH N ONLY IF YOU WISH
          ;TO BACKTRACK TO THE PREVIOUS QUESTION
          ;AND BRANCH FURTHER BACK FROM THERE.
          ;IF THE SIMULATOR ABORTS,
          ;MANUALLY TRIGGER THE SYSTEM
          ;AND ANSWER THE SUBSEQUENT "ABORT?"
          ;WITH Y TO LOOP BACK TO RECOVER
.1010:     .ASK TRIG NEAR TRIGGER TIME
          .IFF TRIG   .GOTO 910
          RUN DATA
          .IFT DIALOG; DATA SHOT TAKEN

```

APPENDIX C

.ASK ABORT ABORT  
.IFT ABORT .GOTO 900

RUN SCALE  
RUN RAWCPY  
RUN MARK  
RUN ERASE  
RUN TABIN  
RUN SCACPY  
RUN TABLES

RUN ERASE  
;ONCE USER COMMENTS HAVE BEEN FILED  
;(SEE XMTR LOG FOR NEXT SHOT) ,  
;RUN THE TASK USRCOM.  
;THE SHOT FILE WILL THEN BE COMPLETE  
;AND SHOULD BE TRANSFERRED  
;TO THE DAILY DATA TAPE WITH FLX.  
;AS AN AID, RUNNING THE TASK PRTNAM  
;WILL GIVE THE ACTUAL NAME  
;OF THE PRESENT SHOT FILE.

# APPENDIX D.--SYMDEF MACRO ENTRIES

SYMDEF is a macro that is called by virtually every program in the SKYNET Applications Software Package (SASP). Its function is to define SASP parameters.

<u>Entry</u>	<u>Assigned value</u>	<u>Definition</u>
ADRCB1	164000	UNIBUS address, Tektronix R7912 digitizer controller
ADRCB2	165000	UNIBUS address, digitizer controller card No. 2
ADRCB3	166000	UNIBUS address, digitizer controller card No. 3
ADRCB4	167000	UNIBUS address, digitizer controller card No. 4
ARMBIT	000020	Digitizer status bit for armed sweep
ATNDLY	000010	Offset for device delay in ATNCAL PSECT (program section)
ATNERR	000004	Bit for transmitter attenuator setting verification error
ATNGN	000004	Offset for device gain in ATNCAL PSECT
ATNID	000000	Offset for device identification (ID) in ATNCAL PSECT
ATNMDL	000014	Offset for device model in ATNCAL PSECT
ATNSER	000024	Offset for device manufacturer's serial number in ATNCAL
BADBIT	010000	Bit to indicate bad (incorrect) digitizer data word; set by software
BLMERR	001000	Bit for digitizer bloom error
BLMFAC	000004	Bloom multiplication factor for average to center of trace (ATC) routines
BLMMAX	000017	Maximum allowed bloom
BLMMIN	000002	Minimum allowed bloom
BZYBIT	020000	Digitizer status bit for busy condition

# APPENDIX D

<u>Entry</u>	<u>Assigned value</u>	<u>Definition</u>
BZYERR	000020	Bit for digitizer busy error
CALERR	000010	Channel error bit for bad (below-standard) calibration trace
CHNATN	000036	Offset for transmitter attenuator setting in channel entry in SYSCON (COMMON area image of SYSCON.DAT file)
CHNAT1	000010	Offset for attenuator ID in channel entry in SYSCON
CHNAT2	000030	Offset for device ID in channel entry in SYSCON
CHNBIS	000020	Offset for x-ray detector (XRD) bias in channel entry in SYSCON
CHNDGS	000052	Offset for list of digitizers associated with channel in channel's entry in SYSCON
CHNDLG	000060	Offset for data link gain in channel entry in SYSCON
CHNDLY	000070	Offset for overall channel delay in channel entry in SYSCON
CHNFLG	000000	Offset for channel error flag in channel entry in SYSCON
CHNNDG	000050	Offset for number of digitizers associated with channel entry in SYSCON
CHNPRA	000044	Offset for postreceiver gain value in channel entry in SYSCON
CHNPRB	000014	Offset for probe/sensor ID in channel entry in SYSCON
CHNPRD	000040	Offset for postreceiver cable delay in channel entry in SYSCON
CHNPXD	000024	Offset for pretransmitter cable delay in channel entry in SYSCON
CHNSFR	000064	Offset for scale factor in channel entry in SYSCON
CHNTPN	000002	Offset for test point name in channel entry in SYSCON
CHNXMT	000034	Offset for transmitter address (or number) in channel entry in SYSCON

<u>Entry</u>	<u>Assigned value</u>	<u>Definition</u>
CH1BIT	000100	Digitizer status bit for channel No. 1 selected (for dual-channel plug-in unit in digitizer)
CLKBGN	000031	Command to start programmable clock for calorimeter data-taking mode
CLRBIT	000216	Digitizer status bits for complement of initialized (cleared) condition
CMDALL	002006	Digitizer command to enable data taking on both channels at dual vertical plug-in
CMDARM	002040	Digitizer command to arm single sweep
CMDCH1	002013	Digitizer command to select channel No. 1 of dual input to digitizer
CMDCH2	002003	Digitizer command to select channel No. 2 of dual input to digitizer
CMDCLR	004000	Digitizer command to initialize registers
CMDDGT	002015	Digitizer command to switch to digital mode
CMDDTA	002100	Digitizer command to start sending data
CMDERR	000040	Bit for digitizer command execution error
CMDGRT	002014	Digitizer command to enable graticule
CMDKRO	002002	Digitizer command to enable sending knob readout information immediately following data
CMDLCK	002017	Digitizer command to lock memory (write-protect)
CMDNGT	002004	Digitizer command to disable graticule
CMDNKO	002012	Digitizer command to inhibit sending of knob readout information
CMDREM	002016	Digitizer command to enable channel selection (dual input only)
CMDSTR	002020	Digitizer command to request status
CMDTV	002005	Digitizer command to switch to television mode

# APPENDIX D

<u>Entry</u>	<u>Assigned value</u>	<u>Definition</u>
CMDUNL	002007	Digitizer command to unlock memory (remove write-protect)
CNTBZY	000036	Number of seconds for time-out on busy status for digitizer
CNTCMD	000036	Number of seconds for time-out on command for digitizer
CNTXMT	000036	Number of seconds for transmitter warm-up clock
CRGCLK	172544	UNIBUS address of count register of programmable clock
CSBCLK	172542	UNIBUS address of count set buffer of programmable clock
CSRCLK	172540	UNIBUS address of command and status register of programmable clock
DGTBIT	000200	Digitizer status bit for digital mode
DIGADR	000000	Offset for digitizer UNIBUS address in digitizer's entry in SYSCON
DIGBAS	000014	Offset for requested baseline trace position in digitizer entry in SYSCON
DIGBLM	000012	Offset for bloom in digitizer entry in SYSCON
DIGCHN	000002	Offset for associated channel number in digitizer's entry in SYSCON
DIGFLG	000004	Offset for digitizer error flag in digitizer's entry in SYSCON
DIGHRZ	000006	Offset for requested horizontal sweep rate setting in digitizer's entry in SYSCON
DIGVRT	000010	Offset for requested vertical sensitivity setting in digitizer's entry in SYSCON
DISERR	000002	Bit for digitizer or channel disabled error
DTAMSK	177000	Mask for extraneous bits in digitizer data word
GRTBIT	000010	Digitizer status bit for graticule enabled

## APPENDIX D

<u>Entry</u>	<u>Assigned value</u>	<u>Definition</u>
HDRASC	000014	Offset for American Standard Code for Information Interchange shot number in SYSCON header
HDRCLE	000144	Offset for calorimeter error flag in SYSCON header
HDRCL1	000130	Offset for channel A calorimeter result in SYSCON header
HDRCL2	000134	Offset for channel B calorimeter result in SYSCON header
HDRCOM	000070	Offset for simulator comment in SYSCON header
HDRCON	000046	Offset for pulser configuration code in SYSCON header
HDRDAT	000000	Offset for trigger date in SYSCON header
HDRIDP	000064	Offset for current diode peak value in SYSCON header
HDRLEV	000050	Offset for pulser output level in SYSCON header
HDRNBK	000126	Offset for number of blocks in MAIN.DAT or shot file in SYSCON header
HDRNCB	000024	Offset for number of digitizer controller boards in system in SYSCON header
HDRNCH	000026	Offset for number of channels in system in SYSCON header
HDRNDG	000022	Offset for number of digitizers in system in SYSCON header
HDRNS1	000140	Offset for number of samples taken on calorimeter channel A in SYSCON header
HDRNS2	000142	Offset for number of samples taken on calorimeter channel B in SYSCON header
HDRSCV	000050	Offset for simulator charging voltage in SYSCON header
HDRSHN	000020	Offset for integer shot number in SYSCON header
HDRSIM	000044	Offset for simulator code in SYSCON header
HDRTIM	000006	Offset for trigger time in SYSCON header

APPENDIX D

<u>Entry</u>	<u>Assigned value</u>	<u>Definition</u>
HDRVAC	000054	Offset for tank vacuum in SYSCON header
HDRVDP	000060	Offset for voltage diode peak value in SYSCON header
HDRWIR	000046	Offset for wire type code in SYSCON header
HZMSK	001000	Bit set to indicate horizontal word in raw digitizer data
INTBIT	000100	Enable bit for digitizer interrupts
KROERR	000400	Bit for digitizer knob readout error
KROMSK	177600	Mask for extraneous bits in digitizer knob readout word
LKDBIT	000004	Digitizer status bit for locked memory (write-protect)
MISMAX	000012	Maximum number of mistakes in digitizer read verification allowed before another reading is requested
NEGMSK	100000	Bit set to indicate knob readout word in digitizer data
NOBERR	020000	Bit set for no digitizer baseline trace error
NOCERR	040000	Bit set for no digitizer calibration trace error
NODERR	100000	Bit set for no digitizer data trace error
NPTMIN	000144	Minimum number of points allowed for correct bloom calculation
NRDBLK	000010	Number of disc blocks required for raw digitizer data trace storage
NSCBLK	000005	Number of disc blocks required for SYSCON
OFFERR	000001	Bit set for digitizer or transmitter off-line error
OVFERR	000004	Bit set for digitizer data overflow error
PRBBLN	000026	Offset for balun number in PRBCAL PSECT
PRBCON	000054	Offset for connector in PRBCAL PSECT

<u>Entry</u>	<u>Assigned value</u>	<u>Definition</u>
PRBDLY	000010	Offset for probe delay in PRBCAL PSECT
PRBFLG	000024	Offset for integration flag in PRBCAL PSECT
PRBID	000000	Offset for probe sensor ID in PRBCAL PSECT
PRBMDL	000034	Offset for model in PRBCAL PSECT
PRBMT	000030	Offset for measurement type in PRBCAL PSECT
PRBPAR	000064	Offset for parameter in PRBCAL PSECT
PRBSER	000044	Offset for serial number in PRBCAL PSECT
PRBTF	000004	Offset for transfer function in PRBCAL PSECT
PRBUNT	000014	Offset for units in PRBCAL PSECT
PRICLK	000300	Priority for clock interrupts
PRIDIG	000300	Priority for digitizer interrupts
PRITRG	000300	Priority for trigger interrupts
RAWMAX	003720	Maximum number of raw digitizer data words allowed
RAWMIN	001750	Minimum number of raw digitizer data words allowed
REMBIT	000002	Digitizer status bit for enabled channel selection
SGLBIT	000040	Digitizer status bit for single sweep mode
SHFERR	002000	Bit set for digitizer baseline trace positioning error
SHFMAX	000015	Maximum allowed number of points between digitizer baseline trace position and requested position
SIZCHN	000100	Size in bytes of channel entry in SYSCON
SIZDIG	000030	Size in bytes of digitizer entry in SYSCON
SIZHDR	000240	Size in bytes of SYSCON header
TANMAX	000080	Maximum allowed absolute value of slope between two points of ATC of digitizer trace
TRGERR	000100	Bit set for digitizer failure-to-trigger error

## APPENDIX D

<u>Entry</u>	<u>Assigned value</u>	<u>Definition</u>
TRYMAX	000012	Maximum allowed number of attempts at verified read of digitizer
UNFERR	000010	Bit set for digitizer data underflow error
UNTMSK	177707	Mask for digitizer unit number in controller status register (upon interrupt) of digitizer controller card
VBNBAS	000030	Virtual block number for beginning of baseline trace storage in MAIN.DAT or shot file
VBNCAL	000032	Virtual block number for beginning of calibration trace storage in MAIN.DAT or shot file
VBNPAR	000040	Virtual block number for beginning of parameter storage in MAIN.DAT or shot file
VBNRAW	000034	Virtual block number for beginning of raw data trace storage in MAIN.DAT or SHOT.DAT file
VBNSCA	000036	Virtual block number for beginning of scaled data trace storage in MAIN.DAT or shot file
VBNUSR	000042	Virtual block number for beginning of user comment storage in MAIN.DAT or shot file
VECCB1	000374	Vector address for digitizer controller card No. 1
VECCB2	000370	Vector address for digitizer controller card No. 2
VECCB3	000364	Vector address for digitizer controller card No. 3
VECCB4	000360	Vector address for digitizer controller card No. 4
VECCLK	000104	Vector address for programmable clock
VECTRG	000334	Vector address for trigger
VFYERR	000200	Bit set for digitizer read verification error
XCHCAL	000014	Offset for calibration generator output level in XCHCAL PSECT
XCHDLY	000010	Offset for delay in XCHCAL PSECT
XCHFBR	000002	Offset for fiber number in XCHCAL PSECT

<u>Entry</u>	<u>Assigned value</u>	<u>Definition</u>
XCHREC	000004	Offset for receiver number in XCHCAL PSECT
XCHXMT	000000	Offset for transmitter address number in XCHCAL PSECT
XHDBLN	000030	Offset for balun insertion loss table in XCHHDR PSECT
XHDDR1	000010	Offset for first part of trigger path delay ( $\Delta t_{\text{trig}_1}$ ) in XCHHDR PSECT
XHDDR2	000014	Offset for second part of trigger path delay ( $\Delta t_{\text{trig}_2}$ ) in XCHHDR PSECT
XHDDTF	000020	Offset for fiducial path delay ( $\Delta t_{\text{fidu}}$ ) in XCHHDR PSECT
XHDDTV	000024	Offset for vertical plug-in delay ( $\Delta t_{\text{vert}}$ ) in XCHHDR PSECT
XHDDT1	000004	Offset for source to XRD transit time in XCHHDR PSECT
XHDNCH	000000	Offset for number of channels entered in XCHCAL.DAT file in XCHHDR PSECT

# APPENDIX E.--LIBRARIES: OBJECT AND MACRO

To facilitate program development, maintenance and understanding, similar subroutines and macros are grouped together into object or macro libraries resident on the system disc.

<u>Name</u>	<u>Contents</u>	<u>Purpose</u>
ERROR (object)	CHNERR DCDBIT DIGERR DISABL DSKERR ERRNUM IASCII LSTERR .RSTR .SAVE	List and flag digitizer and channel errors
HANDLR (object)	BUSY CHKTRG CLRCBD CMDDIG FLGPTS INTABL INTDIG RDIG SETVEC	Digitizer handling routines; many reference error reporting routines in library ERROR
INPUT (object)	ANSWER ASDCD FIELD FLDCD FLOTE FXDCD IFICKS INPUT	Terminal input routines
NORMAL (object)	BASATC CALATC DATATC ELIMPT EXTRAP INTRP ISBASE QKELIM QKLINTP SUBASE	Normalization routines

# APPENDIX E

<u>Name</u>	<u>Contents</u>	<u>Purpose</u>
SKYNET (macro)	RDSECT SCSECT SYMDEF	General, system-wide macros
RRPMAC (macro)	ADR CALL  COM  CMPD  FIX  FLOAT  FPA  FPC  FPE  FPS  FSTK  FSTO  HDR  MIN  MOVF  MOVP  SCRPT  SET  ST  STE  STZ	General purpose macros

# DISTRIBUTION

DEFENSE DOCUMENTATION CENTER  
CAMERON STATION, BUILDING 5  
ALEXANDRIA, VA 22314  
ATTN DDC-TCA (12 COPIES)

COMMANDER  
USA RSCH & STD GP (EUR)  
BOX 65  
FPO NEW YORK 09510  
ATTN LTC JAMES M. KENNEDY, JR.  
CHIEF, PHYSICS & MATH BRANCH

COMMANDER  
US ARMY MATERIEL DEVELOPMENT  
& READINESS COMMAND  
5001 EISENHOWER AVENUE  
ALEXANDRIA, VA 22333  
ATTN DRXAM-TL, HQ TECH LIBRARY

COMMANDER  
US ARMY ARMAMENT MATERIEL  
READINESS COMMAND  
ROCK ISLAND ARSENAL  
ROCK ISLAND, IL 61201  
ATTN DRSAR-ASF, FUZE &  
MUNITIONS SPT DIV

COMMANDER  
USA MISSILE & MUNITIONS  
CENTER & SCHOOL  
REDSTONE ARSENAL, AL 35809  
ATTN ATSK-CTD-F

DIRECTOR  
DEFENSE COMMUNICATIONS AGENCY  
DEPT OF DEFENSE  
WASHINGTON, DC 20305  
ATTN PLANS & PROGRAMS DIRECTORATE  
ATTN OFC OF SYS ENGINEERING  
ATTN TECH LIBRARY  
ATTN NMR

DIRECTOR  
DEFENSE NUCLEAR AGENCY  
WASHINGTON, DC 20305  
ATTN PETER HAAS, DEP DIR,  
SCIENTIFIC TECHNOLOGY  
ATTN RADIATION DIRECTORATE  
ATTN ELECTRONIC VULNERABILITY  
DIV (RAEV)  
ATTN ANALYSIS & PROGRAMS DIRECTORATE  
ATTN TECHNICAL LIBRARY (APTL)  
ATTN VULNERABILITY DIRECTORATE  
ATTN INFORMATION SYSTEMS DIV (VLIS)  
ATTN DDST  
ATTN TITL TECH LIBRARY  
ATTN TISI ARCHIVES

DIRECTOR OF DEFENSE RESEARCH  
& ENGINEERING  
WASHINGTON, DC 20301  
ATTN DEP DIR (TEST & EVALUATION)  
ATTN ASST DIR (ENG TECHNOLOGY)

COMMANDER  
FIELD COMMAND  
DEFENSE NUCLEAR AGENCY  
KIRTLAND AFB, NM 87115  
ATTN FCSD-A4, TECH REF BR  
ATTN D. SHAMBLIN

COMMANDER  
AF WEAPONS LAB, AFSC  
KIRTLAND AFB, NM 87117  
ATTN SE, NUCLEAR SYS DIV  
ATTN EL, ELECTRONICS DIV  
ATTN EL, ELECTRONICS DIV  
ATTN ELC, COMM & SATELLITE BR  
ATTN LR, LASER DEV DIV  
ATTN LRL, LASER DEVICE BR

DIRECTOR  
DEFENSE ADVANCED RESEARCH PROJECT AGENCY  
ARCHITECT BUILDING  
1400 WILSON BLVD  
ARLINGTON, VA 22209  
ATTN NMR

CHIEF  
LIVERMORE DIVISION FIELD COMMAND, DNA  
LAWRENCE LIVERMORE LABORATORY  
P.O. BOX 808  
LIVERMORE, CA 94550  
ATTN FCPRL

COMMANDER  
REDSTONE SCIENTIFIC INFORMATION CTR  
US ARMY MISSILE COMMAND  
REDSTONE ARSENAL, AL 35809  
ATTN CHIEF, DOCUMENTS

CHIEF  
US ARMY COMMUNICATIONS SYS AGENCY  
FORT MONMOUTH, NJ 07703  
ATTN SCCM-AD-SV, LIBRARY

COMMANDER  
US ARMY ELECTRONICS COMMAND  
FORT MONMOUTH, NJ 07703  
ATTN DRSEL

COMMANDER  
US ARMY FOREIGN SCIENCE & TECH CTR  
220 7TH STREET, NE  
CHARLOTTESVILLE, VA 22901  
ATTN DRXST-ISI

DISTRIBURION (Cont'd)

DIRECTOR  
NAVAL RESEARCH LABORATORY  
WASHINGTON, DC 20375  
ATTN CODE 7756, J. DAVIS

COMMANDER  
NAVAL SURFACE WEAPONS CENTER  
WHITE OAK, SILVER SPRING, MD 20910  
ATTN CODE WA501, NAVY NUC PRGMS OFF

AF WEAPONS LABORATORY, AFSC  
KIRTLAND AFB, NM 87117  
ATTN SUL

AF WEAPONS LABORATORY, AFSC  
KIRTLAND AFB, NM 87117  
ATTN NTS

AF WEAPONS LABORATORY, AFSC  
KIRTLAND AFB, NM 87117  
ATTN DYC

COMMANDER  
ROME AIR DEVELOPMENT CENTER, AFSC  
HANS COM AFB, MA 01731  
ATTN EDWARD A. BURKE

SAMSO/DY  
POST OFFICE BOX 92960  
WORLDWAY POSTAL CENTER  
LOS ANGELES, CA 90009  
ATTN DYS (TECHNOLOGY)

SAMSO/MN  
NORTON AFB, CA 92409  
ATTN MNNG (MINUTEMAN)  
ATTN MNNH (MINUTEMAN)

SAMSO/SK  
POST OFFICE BOX 92960  
WORLDWAY POSTAL CENTER  
LOS ANGELES, CA 90009  
ATTN SKF (SPACE COMM SYSTEMS)

SAMSO/XR  
POST OFFICE BOX 92960  
WORLDWAY POSTAL CENTER  
LOS ANGELES, CA 90009  
ATTN XRS (DEVELOPMENT PLANS)

UNIVERSITY OF CALIFORNIA  
LAWRENCE LIVERMORE LABORATORY  
P.O. BOX 808  
LIVERMORE, CA 94550  
ATTN TECH INFO DEPT L-3

NASA  
LEWIS RESEARCH CENTER  
21000 BROOKPARK ROAD  
CLEVELAND, OH 44135  
ATTN LIBRARY (NO R/D)

AEROSPACE CORPORATION  
P.O. BOX 92957  
LOS ANGELES, CA 90009

CALIFORNIA AT SAN DIEGO, UNIVERSITY OF  
IPAPS, B-019  
LA JOLLA, CA 92093  
ATTN SHERMAN DE FOREST

COMPUTER SCIENCES CORPORATION  
201 LA VETA DRIVE N.E.  
ALBUQUERQUE, NM 87108  
ATTN ALVIN T. SCHIFF

DEPLOMB, DR. EUGENE P.  
2422 DEER PARK DRIVE  
SAN DIEGO, CA 92110  
ATTN EUGENE P. DEPLOMB

DIKEWOOD INDUSTRIES, INC  
1009 BRADBURY DRIVE, S.E.  
ALBUQUERQUE, NM 87106  
ATTN TECH LIB

EG&G, INC.  
ALBUQUERQUE DIVISION  
P.O. BOX 10218  
ALBUQUERQUE, NM 87114  
ATTN TECHNICAL LIBRARY

FORD AEROSPACE & COMMUNICATIONS CORP  
3939 FABIAN WAY  
PALO ALTO, CA 94303  
ATTN DONALD R. MCMORROW, MS G30

GENERAL ELECTRIC COMPANY  
SPACE DIVISION  
VALLEY FORGE SPACE CENTER  
GODDARD BLVD KING OF PRUSSIA  
P.O. BOX 8555  
PHILADELPHIA, PA 19101  
ATTN JOSEPH C. PEDEN VFSC, RM 4230M

GENERAL ELECTRIC COMPANY  
TEMPO-CENTER FOR ADVANCED STUDIES  
816 STATE STREET (P.O. DRAWER QQ)  
SANTA BARBARA, CA 93102  
ATTN DASIAC

HUGHES AIRCRAFT COMPANY  
CENTINELA AND TEALE  
CULVER CITY, CA 90230  
ATTN TECH LIB

IRT CORPORATION  
P.O. BOX 81087  
SAN DIEGO, CA 92138  
ATTN TECHNICAL LIBRARY

DISTRIBUTION (Cont'd)

JAYCOR  
1401 CAMINO DEL MAR  
DEL MAR, CA 92014  
ATTN LIBRARY

KAMAN SCIENCES CORPORATION  
P.O. BOX 7463  
COLORADO SPRINGS, CO 80933  
ATTN LIBRARY

MISSION RESEARCH CORPORATION  
735 STATE STREET  
SANTA BARBARA, CA 93101  
ATTN CONRAD L. LONGMIRE

MISSION RESEARCH CORPORATION-SAN DIEGO  
P.O. BOX 1209  
LA JOLLA, CA 92038  
ATTN LIBRARY

R & D ASSOCIATES  
P.O. BOX 9695  
MARINA DEL REY, CA 90291  
ATTN TECHNICAL LIBRARY

SCIENCE APPLICATIONS, INCORPORATED  
8400 WESTPARK DRIVE  
MCLEAN, VA 22101  
ATTN WILLIAM L. CHADSEY

SIMULATION PHYSICS, INC.  
41 "B" STREET  
BURLINGTON, MA 01803  
ATTN ROGER G. LITTLE

STANFORD RESEARCH INSTITUTE  
333 RAVENSWOOD AVENUE  
MENLO PARK, CA 94025  
ATTN LIBRARY

SYSTEMS, SCIENCE AND SOFTWARE, INC.  
P.O. BOX 1620  
LA JOLLA, CA 92038  
ATTN TECHNICAL LIBRARY

TRW DEFENSE & SPACE SYS GROUP  
ONE SPACE PARK  
REDONDO BEACH, CA 90278  
ATTN ROBERT M. WEBB, RL-2410

HARRY DIAMOND LABORATORIES  
ATTN RAMSDEN, JOHN J., LTC, COMMANDER/  
FLYER, I.N./LANDIS, P.E./  
SOMMER, H./OSWALD, R. B.  
ATTN CARTER, W.W., DR., TECHNICAL  
DIRECTOR/MARCUS, S.M.  
ATTN WISEMAN, ROBERT S., DP.,  
DRDEL-CT  
ATTN KIMMEL, S., PAO  
ATTN CHIEF, 0021  
ATTN CHIEF, 0022  
ATTN CHIEF, LAB 100  
ATTN CHIEF, LAB 200  
ATTN CHIEF, LAB 300  
ATTN CHIEF, LAB 400  
ATTN CHIEF, LAB 500  
ATTN CHIEF, LAB 600  
ATTN CHIEF, DIV 700  
ATTN CHIEF, DIV 800  
ATTN CHIEF, LAB 900  
ATTN CHIEF, LAB 1000  
ATTN RECORD COPY, BR 041  
ATTN HDL LIBRARY (5 COPIES)  
ATTN CHAIRMAN, EDITORIAL COMMITTEE  
ATTN CHIEF, 047  
ATTN TECH REPORTS, 013  
ATTN PATENT LAW BRANCH, 071  
ATTN GIDEP OFFICE, 741  
ATTN LANHAM, C., 0021  
ATTN ROSADO, J., 210  
ATTN BLACKBRUN, J., 230  
ATTN PUTTCAMP, R., 230 (20 COPIES)  
ATTN TOMPKINS, J., 230  
ATTN SCHALLHORN, D., 230  
ATTN PETTY, W., 1050 (2 COPIES)  
ATTN GILBERT, R., 230  
ATTN SCOTT, W., 1050  
ATTN MATTHEWS, H. J., III, 0025  
ATTN CALDWELL, P. A., CHIEF 290  
ATTN LAMB, R. C., 290